Automatic Evaporation

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 Datum

 Unterschrift, Firmenstempel

Translation of the
original instructions

Contents

Introduction

| | |
|------------------------------|---|
| About this document..... | 5 |
| Typographic conventions..... | 5 |
| Copyright protection..... | 5 |

General notes

| | |
|---|---|
| Basic product information | 6 |
| Guidelines applied, product certification | 6 |
| California Residents..... | 6 |
| Residual risk | 6 |
| Intended use..... | 6 |
| Reasonably foreseeable misuse..... | 6 |
| Transportation..... | 7 |
| Storage..... | 7 |
| Acclimatization | 7 |
| Permissible ambient conditions..... | 7 |

Safety

| | |
|--|----|
| General safety information..... | 8 |
| Electrical safety..... | 8 |
| Operational safety | 8 |
| Biohazard | 9 |
| Special hygiene measures for the use of laboratory equipment in food, cosmetics and pharmaceutical production..... | 9 |
| General Measures | 9 |
| Device-specific measures | 10 |
| Other regulations | 10 |

Assembly

| | |
|----------------|----|
| Assembly | 11 |
|----------------|----|

Device description

| | |
|------------------------------|----|
| Principle of operation | 12 |
|------------------------------|----|

| | |
|---|----|
| Mechanical design | 13 |
| Front view control box with control panel | 13 |
| Connections rear side control box..... | 14 |
| Glassware..... | 15 |
| Control panel..... | 16 |
| Structure of the user interface..... | 17 |
| Control and display elements..... | 18 |
| System time/Timer settings..... | 19 |
| Color coding..... | 20 |
| Symbols | 20 |
| Input keys | 21 |
| User input..... | 21 |
| Commissioning | |
| Commissioning..... | 22 |
| Switch the device on/off | 22 |
| Calibrate the flask level sensor | 22 |
| Initial filling of the evaporation flask | 25 |
| Device configuration..... | 26 |
| Hauptmenü | 26 |
| Applications menu | 28 |
| Settings menu | 29 |
| Screen lock | 30 |
| Display brightness..... | 31 |
| Calibrate filling level | 31 |
| Factory settings | 31 |
| Select the user language | 31 |
| Select formats | 32 |
| Automatic drainage of residue | 33 |
| Select evaporator type/glassware..... | 33 |
| Data Recording | 34 |
| System check | 36 |
| Test valves and signaling..... | 37 |
| Perform overpressure test..... | 38 |
| Test the safety switch | 40 |

Operation

| | |
|------------------------------------|----|
| Select operating mode | 41 |
| Sensor mode | 41 |
| Parameter sensor mode..... | 45 |
| Parameters Hei-VAP | 48 |
| Data logging | 48 |
| Recipes..... | 49 |
| Recipe list..... | 50 |
| Time mode | 51 |
| Time mode parameters..... | 53 |
| Rinse mode | 56 |
| Rinse mode parameters | 58 |
| Manual mode | 61 |
| View error list | 63 |
| Restore factory settings..... | 63 |
| Access system information..... | 64 |
| Firmware and software update | 64 |

Troubleshooting

| | |
|----------------------|----|
| Troubleshooting..... | 65 |
|----------------------|----|

Attachments

| | |
|---|----|
| Technical Specifications | 67 |
| Scope of delivery..... | 68 |
| Device service..... | 69 |
| General cleaning instructions..... | 69 |
| Repairs..... | 69 |
| Maintenance..... | 69 |
| Dismantling..... | 70 |
| Disposal | 70 |
| Contact details Germany – Austria – Switzerland | 71 |
| Warranty Statement | 71 |
| Declaration of no objection | 72 |
| RoHS Declaration of Conformity..... | 73 |
| China RoHS Certification..... | 74 |

About this document

This operating manual describes the features and operation of automatic modules type **Hei-VOLUME Distimatic Pro**. The operating instructions manual is an integral part of the delivery!

Typographic conventions

Standardized symbols, highlighting elements, and signal words are used in this document to identify warnings, cautions, important information, and special text contents.

| Symbol | Signal word / explanatory note |
|---|---|
|  | <p>Warning symbols in combination with a signal word indicate dangers:</p> <p>DANGER Indicates a hazardous situation which, if not avoided, will result in death or serious injury.</p> <p>WARNING Indicates a possible hazardous situation which, if not avoided, could result in death or serious injury.</p> <p>CAUTION Indicates a hazardous situation which, if not avoided, could result in minor or moderate injury or material/environmental damage.</p> |
|  | <p>Prohibition signs indicate actions or situations that are to be omitted or avoided. Failure to comply may result in personal injury and / or damage to property.</p> |
|  | <p>Mandatory signs are used to indicate important information regarding the product handling. This information is used to ensure operational safety and to maintain the value of the product.</p> |
| [GUI] | <p>Parameter designations, display texts, and device labels are highlighted in text and tables in a typographic manner to facilitate the assignment on the device.</p> |
| → | <p>The arrow symbol indicates instructions to be followed in order to ensure the operational safety when handling the product.</p> |

Copyright protection

This publication is protected by copyright and intended for internal use by the purchaser of the product only.

No part of this publication may be transmitted or reproduced in any form, by any means, without the prior written consent of the copyright owner Heidolph Instruments GmbH & Co. KG. Any violation is subject to compensation for damages.

Basic product information

Guidelines applied, product certification

CE Marking

The device complies with the following standards:

- Low Voltage Directive 2006/95/EC
- EMC Directive 2004/108/EC

California Residents

Important information for California residents regarding Prop 65. Please visit www.P65Warnings.ca.gov for more information.

Residual risk

The device was designed and manufactured in accordance with the latest technical standards at the time of development and the recognized safety regulations. During installation and use, as well as during maintenance work, repairs and cleaning, there are nevertheless certain residual risks associated with the device described.

These are identified and described at the appropriate points in this document.

Intended use

The **Hei-VOLUME Distimatic Pro** automatic module, according to the EU Machinery Directive 2006/42/EC is a partly completed machine that has been specially developed for the automatic filling and draining of Heidolph rotary evaporators.

The following processes can be controlled with the automatic module:

- filling educt
- drainage of distillate
- drainage of condensate
- drainage of residue (optional)

Reasonably foreseeable misuse

For use under conditions or for purposes deviating from the intended use, additional measures may become necessary, and/or specific guidelines and safety regulations will have to be observed (see section "Special hygiene measures for the use of laboratory equipment in food, cosmetics and pharmaceutical production" on page 83). Corresponding requirements must be evaluated and observed by the operator in each individual case.

Compliance with and implementation of all relevant guidelines and safety measures for the respective field of application is within the sole responsibility of the operator.

All risks resulting from improper use are solely borne by the operator.

The device may exclusively be operated by authorized and instructed personnel. Training and qualification of the operating personnel as well as ensuring that the device is operated with responsibility are the sole responsibility of the operator!

Transportation

During transport, avoid severe shocks and mechanical stresses that can cause damage to the device.

Keep the original packaging in a dry and protected place for later use.

Storage

Always store the device in its original packaging. To protect against damage and unreasonable material aging, store the device in a dry environment that should be as temperature-stable and dust-free as possible.

Recommended ambient conditions for storage:

- 5 °C – 31 °C up to 80 % rel. humidity
- 32 °C – 40 °C up to 50 % rel. humidity

Acclimatization

After each transport and after storage under critical climatic conditions (e.g. high temperature difference between inside and outside), allow the device to acclimatize at room temperature for a minimum of two hours to prevent possible damage from condensation before putting it into operation at the place of use. If necessary, extend the acclimatization phase if the temperature differences are very high.

Make all supply connections (power supply, tubing) only after the device has been acclimatized!

Permissible ambient conditions

The device is designed for indoor use only. Permissible ambient conditions for operation:

- 5 °C – 31 °C up to 80 % rel. humidity
- 32 °C – 40 °C up to 50 % rel. humidity
- Maximum height above sea level: 2,000 m

When used in corrosive atmospheres, the service life of the device may be reduced depending on the concentration, duration and frequency of exposure.

The device **IS NOT** suitable for outdoor use!

The device **IS NOT** suitable for use in hazardous areas!

General safety information

- Before commissioning and using the device, familiarize yourself with all the safety regulations and guidelines for occupational safety applicable at the place of use and observe them at all times.
- Only operate the device if it is in perfect technical condition. In particular, ensure that there is no visible damage on the device itself and, where necessary, on connected devices or the supply connections.
- If there is missing or misleading information on the device or on occupational safety, contact the responsible safety specialist or our technical service.
- Only use the device in accordance with the regulations for intended use ("Intended use" on page 80).

Electrical safety

- Ensure that the voltage indicated on the rating plate matches the supply voltage of the country in which the device is being used.
- Ensure that the mains socket-outlet is protected by means of a residual-current device (RCD).
- Always use the supplied power supply 3-pin cord provided with the device (phase, neutral, ground).
- Prior to use, check that the device and the power supply cord are free of visible damage.
- Have repairs and/or maintenance work on the device carried out exclusively by an authorized and skilled electrician or by the technical service department of Heidolph Instruments.
- Always switch the device OFF and disconnect it from the power supply before carrying out maintenance work, cleaning, or repairs.

Operational safety

- Do not make any unauthorized changes or modifications to the device!
- Only use genuine spare parts and accessories, or those expressly approved by the manufacturer!
- Rectify malfunctions or faults on the device immediately.
- Shut down the device and disconnect it from the power supply if it is not possible to eliminate the malfunction or rectify the fault immediately.
- Observe all other applicable regulations such as laboratory and workplace guidelines, recognized safety technology rules and special local regulations.
- Observe all relevant general and safety instructions for the connected peripheral devices (observe the supplied documentation!).

Noncompliance will invalidate any warranty against Heidolph Instruments.

The operator is solely liable for all damage resulting from unauthorized changes or modifications to the unit, from the use of unauthorized or non-genuine spare parts and accessories, or from disregarding the safety instructions and hazard warnings or the manufacturer's instructions!

Biohazard

When processing biohazardous substances, take appropriate measures to prevent hazards to persons and the environment, including:

- Instruction of the personnel regarding the necessary safety measures.
- Provision of personal protective equipment (PPE) and instruction of the personnel in its use.
- Marking of the device with the biohazard warning symbol:

The evaluation of corresponding measures such as the marking of a hazardous area, their implementation, and the training of the responsible personnel is the sole responsibility of the operator!

Special hygiene measures for the use of laboratory equipment in food, cosmetics and pharmaceutical production

When laboratory equipment is used in the production processes of the food, cosmetics or pharmaceutical industry, special hygiene measures must be taken by the user to avoid sample contamination and to minimize any risk to humans and the environment as far as possible.

Please observe the following recommendations:

General Measures

- Ensure a clean working and storage environment when handling substances and materials.
- Train all employees in the field of occupational hygiene, document all training measures and check the implementation of all required hygiene measures during operation regularly.
- Use a hygiene control concept such as HACCP (Hazard Analysis and critical Control points). The HACCP comprises the following criteria:
 - Hazard analysis
 - Identification of critical control points
 - Definition of critical limit values
 - Establishment of a system for monitoring and controlling critical hazard control points (CCP)
 - Corrective actions for uncontrollable CCP
 - Establishment of a system to verify the implementation of all HACCP measures
 - Establishment of a system for documenting all associated procedures and protocols

The evaluation of the applicability of the mentioned rules and regulations is within the sole responsibility of the operator!

Device-specific measures

- Regularly clean components that come into contact with the product, such as flasks, seals, tubes, etc. in the autoclave (if available or possible) or chemically (e.g. with ethanol) to sterilize all surfaces.
- Make sure that even products that are intended for single use only are of sufficient purity.
- Do not use open containers.
- Avoid contamination by handling contaminated vessels, apparatus or aids with care.

Contact information

For further information, please contact our after sales service at any time.

Phone: +49-9122-9920-0

Mail: sales@heidolph.de

Other regulations

In addition to the notes and instructions in this document, observe all other applicable regulations such as laboratory and workplace guidelines, hazardous substances ordinances, recognized rules of safety engineering and occupational medicine as well as particular local regulations!

Assembly

If required, Heidolph offers a professional full assembly and commissioning service, including the device and all supplied accessories.

Most of the sensors and valves of the automatic module are plugged into the coded connectors on the back of the control box (sensor cards).

Thanks to the coding of all connectors, all units are automatically detected and assigned in the device software!

All other units must be plugged into the corresponding marked connection sockets and terminals. For a detailed description of all connections, please refer to the assembly instructions supplied.

Faulty or disconnected units are detected and displayed when performing a system test (integrated software function).

The supplied assembly instructions describe all steps for the correct connection of all tubes and cables between the control box, the components supplied (vessels, sensors, valves, glassware) and the evaporation system in detail!

CAUTION: Material damage, loss of production

Due to improper system tubing and wiring as well as unauthorized changes to the system, there is a risk of direct and indirect damage to property!

- For system tubing and wiring, follow all notes and instructions provided in the assembly instructions provided (part number 01-001-009-12-0).
- Do not make any unauthorized changes to the device, system, tubing, nor the existing cable connections.
- If necessary, contact our technical service (see "Contact information Heidolph international" on page 145).

The professional mounting and installation of the device including the supplied accessories is within the sole responsibility of the operator!

Noncompliance will invalidate any warranty against Heidolph Instruments.

The operator is solely liable for all damage resulting from unauthorized changes or modifications to the unit, from the use of unauthorized or non-genuine spare parts and accessories, or from disregarding the safety instructions and hazard warnings or the manufacturer's instructions!

Principle of operation

The **Hei-VOLUME Distimatic Pro** automatic module has been specially developed for use in combination with rotary evaporators. The automatic module enables the automatic filling and draining of the system in continuous operation.

The following processes can be controlled with the **Hei-VOLUME Distimatic Pro** automatic module:

- filling the evaporation flask
- drainage of distillate
- drainage of residue

The following operating modes are available:

- **SENSOR MODE**
Standard mode for automatic filling and draining in continuous operation (see section "Sensor mode" on page 115).
- **TIME MODE**
The time mode is used specifically for processing media that forms sticky residue or that tends to foam (see section "Time mode" on page 125).
- **RINSE MODE**
The rinse mode is used to flush the supply lines and drains with a solvent for cleaning purposes (see section "Rinse mode" on page 130).
- **MANUAL MODE**
In manual mode, the connected rotary evaporator can be operated for processing small quantities (see section "Manual mode" on page 135).

In automatic operation (sensor and time mode), the distillate and – depending on the variant (with or without automatic drainage of residue) – the residue is drained from the collector or the evaporative flask applying overpressure. Distillate and residue can also be drained manually using the function buttons of the HMI. The required overpressure is generated by a compressor integrated into the control box of the automatic module.

The evaporation flask is filled from the educt container using the vacuum in the rotary evaporator system. The required vacuum is generated by an external vacuum pump.

A switching valve is used to switch between the compressor (pressurized operation) and the vacuum pump (vacuum operation).

Mechanical design

This section provides an overview of the mechanical design of the control box and the control panel of the **Hei-VOLUME Distimatic Pro** automatic module. All other components of the delivery vary depending on the product variant and are described in more detail in the respective assembly instructions supplied.

The control box is housed in a sheet steel casing with protection class IP 42 with removable side covers and a removable control panel.

The control box contains the main circuit board, all function modules and the compressor (without illustration) for generating the overpressure to drain the residue and the distillate.

Front view control box with control panel


Connections rear side control box

The following figure shows the rear side of the control box with standard plug-in cards and optional expansion cards. The actual assignment depends on the scope of the order!

The sensors and valves of the system are plugged into the coded connection sockets of the sensor cards. Thanks to the coding of the sockets, all units are automatically detected and assigned in the device software.


In addition, observe the labels of the existing connection sockets and terminals, as well as any kind of information on replacing existing or installing optional expansion cards in the supplied assembly instructions (item no. 01-001-009-12-0) or the supplied documentation for each expansion card!

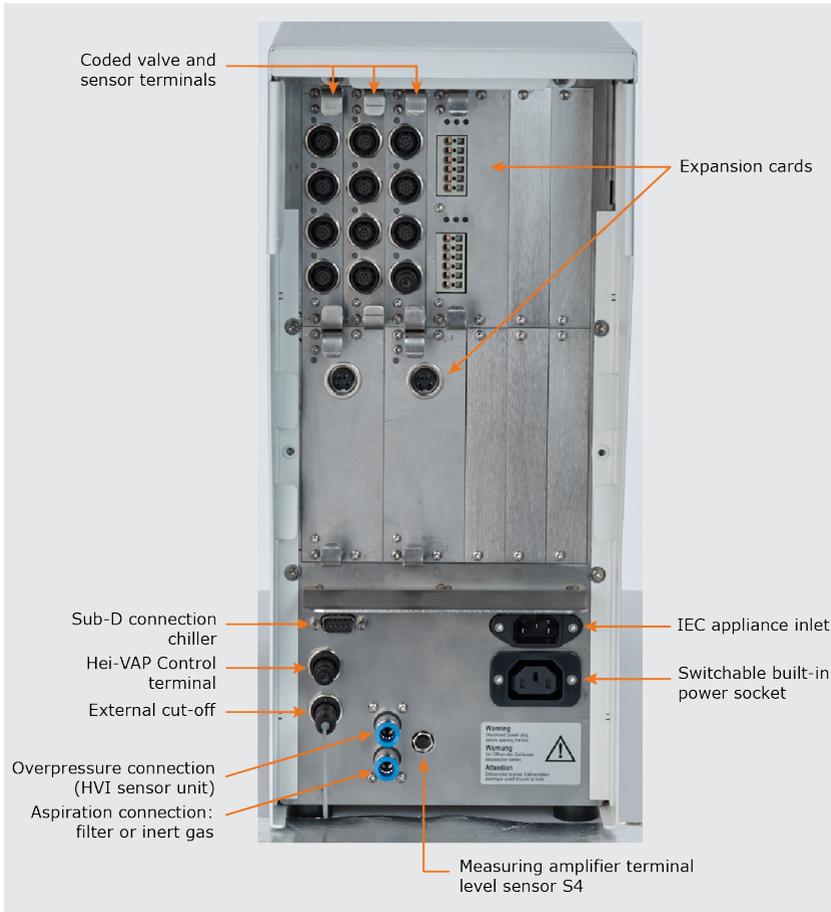Glassware

The use of the automatic module **Hei-VOLUME Distimatic Pro** requires a special set of glassware consisting of a Woulff bottle, a collector vessel and an educt sensor vessel. The required set of glassware is included in the scope of delivery, see section "Scope of delivery" on page 142 and delivery documents.

The following figure shows the maximum configuration with Woulff bottle, collector vessel and educt sensor vessel, including all valves and sensors, completely mounted on the housing of a rotary evaporator:

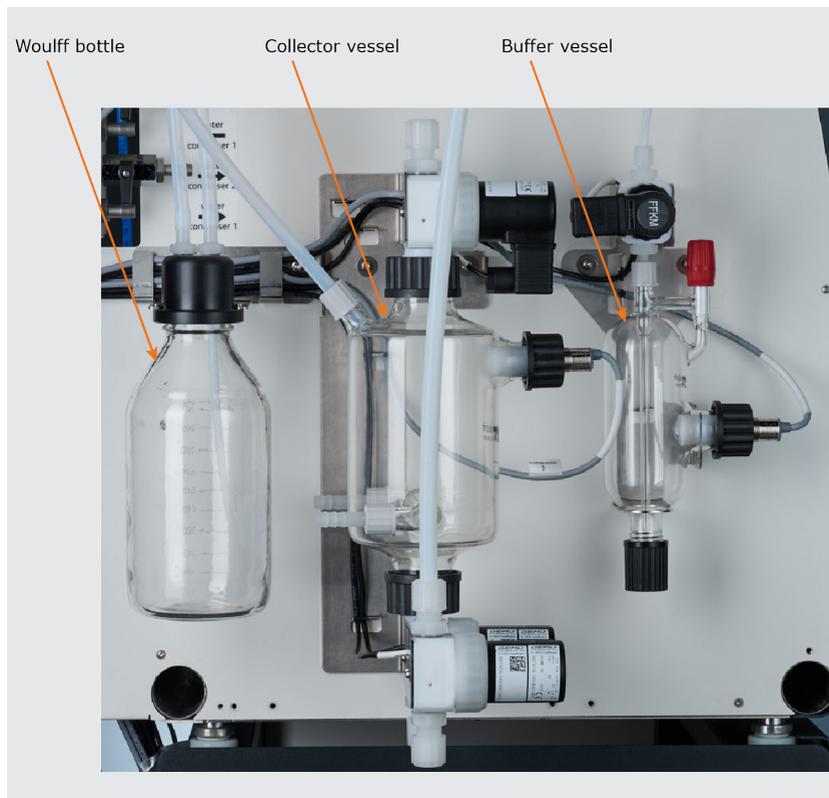When retrofitting an existing system or when installing the set of glassware consisting of Woulff bottle, collector vessel and educt sensor vessel, observe the instructions contained in the assembly instructions delivered (article no. 01-001-009-12-0), as well as the information given in section "Assembly" on page 85 of these operating instructions.

Control panel

The following figure shows the removable control panel of the **Hei-VOLUME Distimatic Pro** automatic module. The control elements of the central display are touch-sensitive and can be operated by finger or with a suitable stylus.

CAUTION: Damage to the control panel

Tip or sharp objects can damage the surface of the display on the control panel!


→ To operate the device, touch the display only with your fingertips or a special stylus with a soft-touch tip.

Structure of the user interface

All device functions are controlled via the integrated control panel and the graphical user interface.

After switching on the device and initializing the device control, the display shows the start window of the last active operating mode (in the example: Sensor mode):

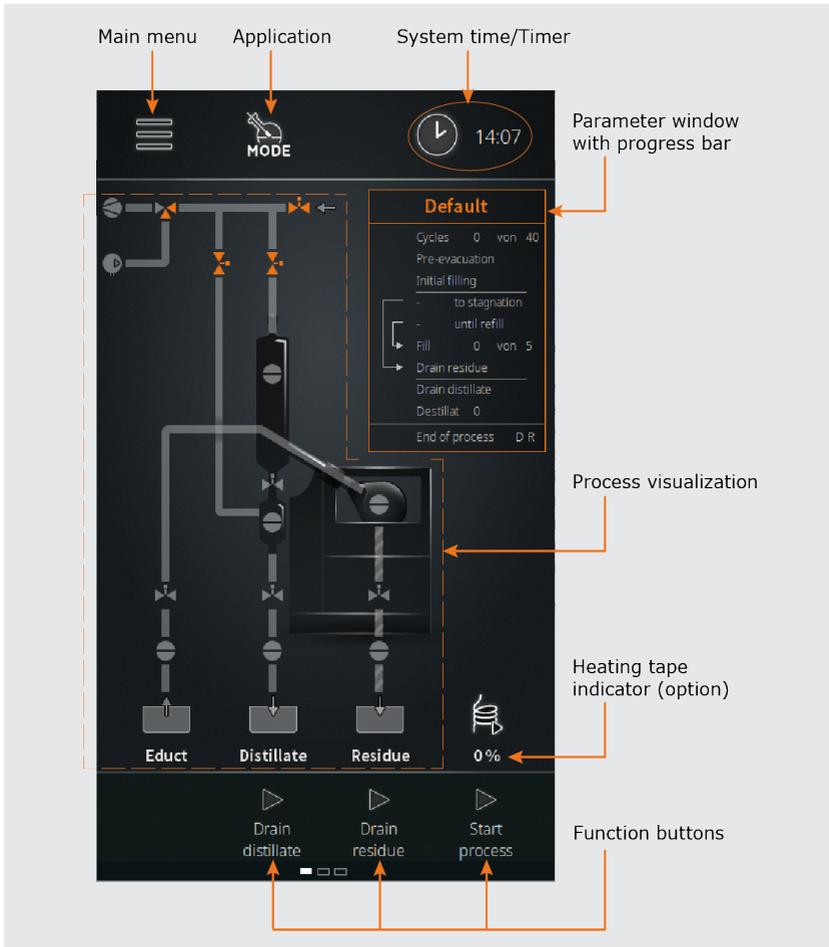From the start window, you can access the main menu, select an operating mode directly, and select the time display, see also "Control and display elements" on page 92.

The function keys in the footer can be used to drain the distillate, to drain the residue, and to start the process directly.

Control and display elements

| Control element | Function |
|---|---|
| Function button Main menu | <p>Touch this function button to open the main menu with the following options:</p> <ul style="list-style-type: none"> ▪ Applications → operating mode ▪ Settings → system parameters ▪ Languages → menu language ▪ Formats → system configuration and display formats ▪ Data logging → data recording ▪ System check → complete system check ▪ Error list → open error list ▪ Info → open system info window |
| Function button Mode | <p>Touch this function button to open the Applications sub-menu directly. Operating modes:</p> <ul style="list-style-type: none"> ▪ Sensor mode ▪ Time mode ▪ Rinse mode ▪ Manual mode |
| Function button System time/Timer | <p>Touch this function button to toggle between the following functions:</p> <ul style="list-style-type: none"> ▪ System time display ▪ Timer ▪ Stopwatch <p>See the following section "System time/Timer settings" on page 93 for detailed information.</p> |
| Drain distillate button | <p>Touch this button to drain the distillate from the collector → only when no process is active!</p> <p>By touching this button, valves V2 and V4 are opened, see section "Principle of operation" on page 86.</p> |
| Drain residue button | <p>Touch this button to drain the residue from the evaporation flask → only when no process is active!</p> <p>By touching this button, valve V3 is opened, see section "Principle of operation" on page 86.</p> |
| Start process button | Touch this button to start the selected process. |
| Default window | The progress of the individual process steps is visualized in this window. |
| Process visualization main display | The current process step is visualized in the main display. |
| Left push-and-rotary control | <ul style="list-style-type: none"> ▪ LED ring off → press briefly to start the process ▪ LED ring on → press long for process stop ▪ The color of the LED ring indicates the process status: green → overpressure, blue → vacuum, white → ventilation & atmospheric pressure |

Right push-and-rotary control

- LED ring lights up orange when activated
 - Function: Increase/decrease value, character forward/backward, confirm value by pressing.
- Press and hold both push-and-rotary controls simultaneously for approximately three seconds to lock the display.

Standby button


Press the standby button to put the device in standby mode. Press the button again to put the device into the previously active operating mode.

Pressing the standby button will stop the current process!

System time/Timer settings

Touch the System time/Timer function button in the header of the start window for about two seconds to open the selection window. In this window, you can set the system time and select the desired display mode. The following options are available:

- Display of the current time (12- or 24-hour display, depending on the selection in the **Settings** menu).
- Stopwatch mode: In this mode, the actual elapsed time from the start of the process is displayed.
- Timer mode: In this mode, the time remaining until the scheduled end of the process is displayed.

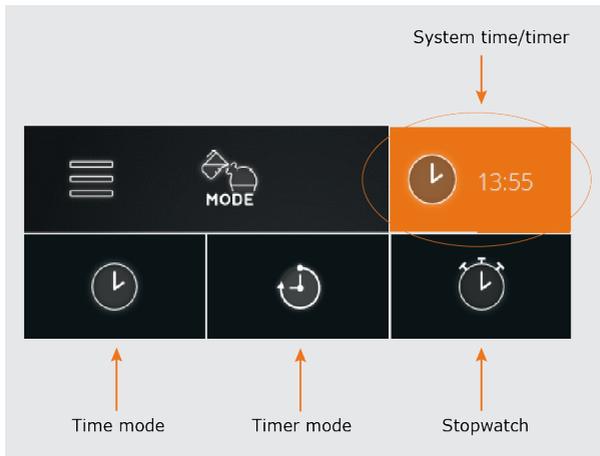This setting is applied in all operating modes and can be changed at any time!

Color coding

In the process visualization area (see section "Structure of the user interface" on page 91), all process lines and function elements are color-coded to illustrate the operating states:

| Color | Meaning |
|------------|---|
| Gray | PASSIVE or off state (e.g. non-active/closed line, non-switched valve, sensor low) |
| Orange | ACTIVE or on state (e.g. active/open line, switched valve, sensor high) |
| White | VENTILATED condition, atmospheric pressure present |
| Light blue | VACUUM condition, vacuum pump active, evacuation |
| Green | OVERPRESSURE condition, compressor on, draining |

Symbols

In the process visualization (see section "Structure of the user interface" on page 91), the following symbols are used to illustrate the process sequences and the integrated function elements (for color coding, see section above):

| Symbol | Meaning |
|---|---------------------|
|  | Compressor |
|  | Vacuum pump |
|  | Pressure mode valve |
|  | 2-way valve |
|  | Sensor |
|  | Heating tape |


Input keys

The following input keys are available in any menu:

| Symbol | Meaning |
|---|---|
|  | Delete button. Press this button to correct a numeric or alphanumeric value entered. |
|  | ok / confirm button. Press this button to confirm and accept a selection or input value. |
|  | Discard button. Press this button to discard a selection or input value and return to the previous view. |

User input

User input for configuration, parameterization and programming is carried out via the numerical and the alphanumeric on-screen keyboard. The on-screen keyboard appears automatically as soon as an input field is touched:

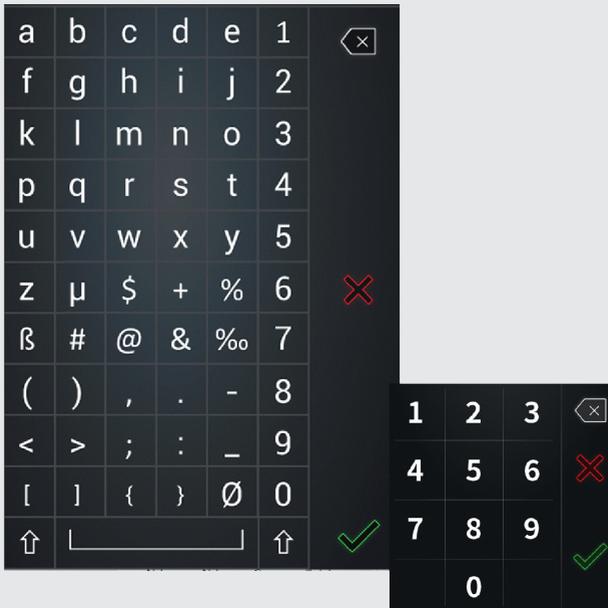Commissioning

Observe the following steps when commissioning the **Hei-VOLUME Distimatic Pro** automatic module.

Switch the device on/off

To switch on and off the automatic module, use the main switch on the front of the control box (see section "Front view control box with control panel" on page 87).

CAUTION: Property damage, loss of production

If the control box is switched on and off incorrectly, there is a risk that peripheral devices will not be detected correctly or will switch in an uncontrolled manner!

- Observe the order given when switching the control box and the connected peripheral devices on and off.
- In particular, before turning off the control box, make sure that all processes have been completed properly.

Observe the following sequence when turning the automatic module and attached peripherals on and off:

- **TURNING ON THE SYSTEM: First, turn on all peripheral devices** so that they are recognized when initializing the automatic module software.
- **TURNING OFF THE SYSTEM: First, turn off all peripherals** to put the entire system in a safe operating condition.

Calibrate the flask level sensor

The flask level sensor (S4) must be calibrated at the first start-up and every time the solvent is changed.

CAUTION: Loss of production

The function of the level sensor in the evaporation flask (S4) may be affected by electromagnetic fields. This can lead to a loss of the operational quality of the entire system.

- During calibration, commissioning and operation, ensure that there are no sources of interference in the immediate vicinity of the evaporation flask to ensure proper operation of the level sensor at all times.

Overflow of the educt is prevented by the emergency-stop sensor S3, see section "Principle of operation" on page 86!


As soon as this sensor responds in the cooler, the control box of the automatic module and the connected peripheral devices are put into an operating safe state and no further medium is fed in.

To calibrate the flask level sensor, follow these steps:

- Switch on the control box and wait approximately 10 minutes for the electronics to warm up to operating temperature. In this way, measurement inaccuracies are avoided during calibration!

Calibration to maximum sensitivity, low-polar solvents:

- Remove the sensor from the evaporation flask.
- Hold the sensor in the open space.
- Make sure there are no objects or bodies in a hemispherical space one (1) meter in front of the sensor head. These would be detected as a signal and falsify the calibration result!
- Start the calibration process in the device software:
 - Open the main menu using the **Main menu** function button in the header of the start window:

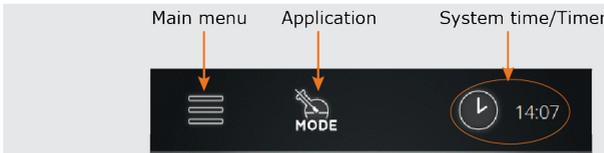
- Open the **Settings** menu and select **Calibrate filling level** (see also "Main menu" on page 100).
- Follow the instructions on the display.

OR

Calibration on high-polar solvents and mixtures:

- Make sure the sensor is in the flask.
- Fill the evaporation flask with at least 500 ml of the product.
- Slide the sensor into the standing piston until the sensor head is just above the fluid level or – for example, in the case of mixtures with water – until the sensor head touches the fluid level slightly.
- Start the calibration process in the device software:
 - Open the main menu using the **Main menu** function button in the header of the start window:

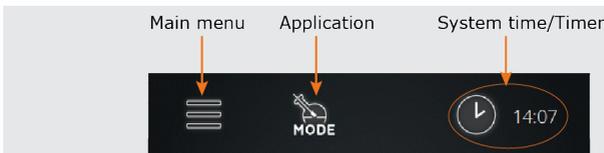- Open the **Settings** menu and select **Calibrate filling level** (see also "Main menu" on page 100).
- Follow the instructions on the display.

PRACTICAL TIP

When processing solvents with a tendency toward boiling retardation (e.g. ethanol), it is recommended to calibrate the sensor slightly less sensitively.

This prevents the boiling delay from being detected as a fluid and that the fluid level in the rotary piston decreases excessively before adding new media.

The sensor should only react to liquids, never to boiling delay or foam!

Calibrate flask level sensor

When calibrating the flask level sensor, the polarity of the medium must be observed! For details see operating instructions, section 'Commissioning - Calibrate flask level sensor'.

Start calibration process with OK!

Calibrate flask level sensor

Calibrating flask level sensor.

Please wait until the indicator changes to green!


- Start the calibration process with .
- The calibration process can be canceled at any time with .
- The progress is visualized by the traffic light display. When the calibration is complete, the green dot on the right lights up.
- Confirm the completion of the calibration process with .
- After the sensor calibration has been completed, the start screen of the last active operating mode is displayed.

Initial filling of the evaporation flask

Regardless of the operating mode selected, the following steps must be followed to fill the evaporation flask for the first time to ensure continuous feed of the educt during operation:


- Turn on the rotary evaporator and wait for the software to initialize (start screen appears on the control panel).
- Switch on the heating bath heater and the recirculation chiller.
- Turn on control box and wait for the software to initialize (start screen appears on the control panel).
- Decrease the system pressure, until a vacuum pressure of 7 mBar is reached (observe the control panel of the rotary evaporator!).
- Start any given process by touching the **Start process** button on the control panel of the Distimatic Pro control box.
- Open the inlet valve on the educt sensor vessel.
- Fill the educt sensor vessel to the bottom edge of the inlet

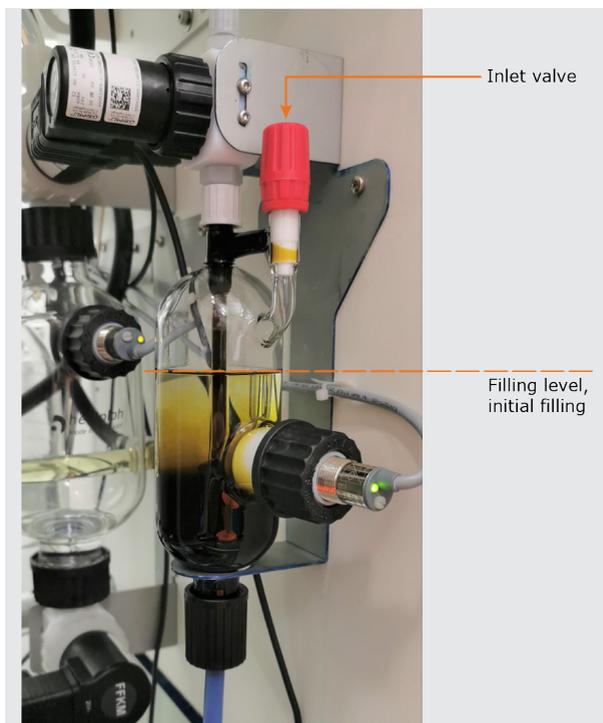
- Stop the process (function button **Start process**, control panel control box Distimatic Pro) to stop the filling.
- If the level in the educt sensor vessel is too high, open the residue valve V3 briefly to lower the level as described.
- Close the inlet valve.

Device configuration

The basic settings for operation are made in the main menu. After switching on the control box and initializing the software, the start page of the last active mode appears. To open the main menu, touch the **Main menu** function button in the header of the display:

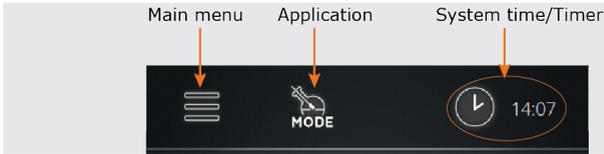Main menu

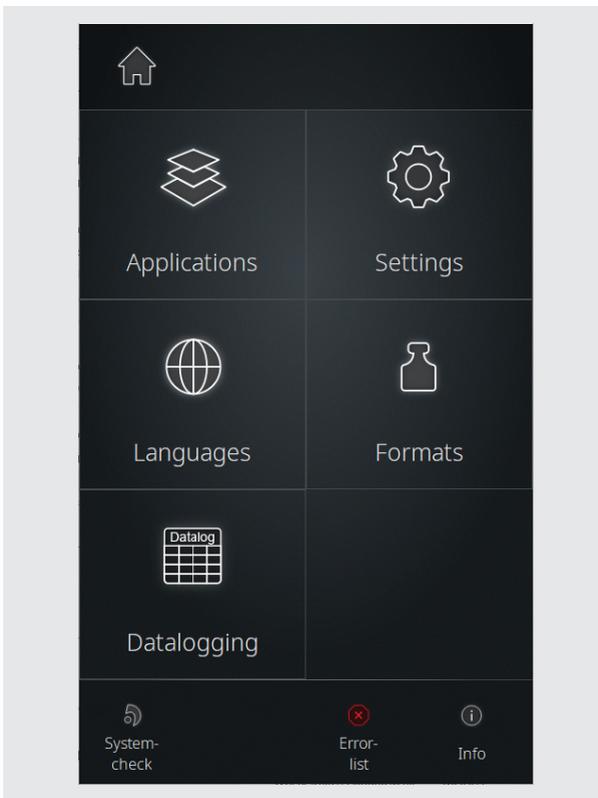
The following options are available from the main menu:

| Button | Function/Meaning |
|---|--|
|  | <p>Home button. Touch this function button to return to the previously opened view (start window).</p> |
| Operating modes | <p>Touch this function button to open the Applications menu and select an operating mode, see section "Applications menu" on page 102.</p> |
| Settings | <p>Touch this function button to open the Settings menu. Various basic settings can be adjusted here and the automatic module can be reset to the factory settings, see section "Settings menu" on page 103.</p> |
| Languages | <p>Touch this function button to open the Languages menu and load a desired user language, see section "Select the user language" on page 105.</p> |
| Formats | <p>Touch this function button to open the Formats menu. Here you can select the desired format for the display of numerical values, time, date and temperature. At the same time, the connected rotary evaporator / glassware is assigned and the operation with or without residue drainage is selected, see section "Select formats" on page 106.</p> |
| Data logging | <p>Touch this function button to start the Data logging function, see section "Data Recording" on page 108.</p> |
| System check | <p>Touch this function button to perform a full system check, see section "" on page 109.</p> |
| Error list | <p>Touch this function button to open the list of all error messages, see section "View error list" on page 137.</p> |
| Info | <p>Touch this function button to open the information window of the system, including firmware and hardware data, see section "Access system information" on page 138.</p> |

Applications menu

The desired operating mode can be selected in the **Applications** menu (**Main menu** → **Applications**).

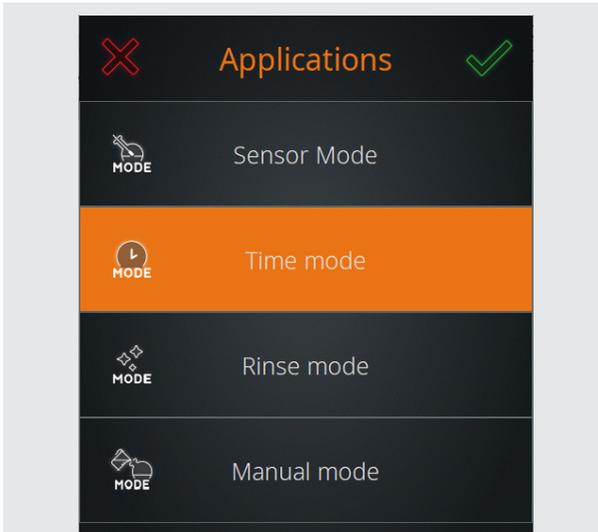To select and activate an operating mode, proceed as follows:


- Complete all running processes before changing the operating mode.
- Touch the entry of the desired operating mode. The background color of the selected entry changes to orange, as shown.
- Press  to confirm and to activate the selected mode.
 - The start window of the activated mode appears on the display.
- With  you can discard the selection and return to the previous view.

The **Applications** menu can also be opened directly by touching the **Mode** function button in the header of the start window.

Settings menu

In the **Settings** menu (**Main menu** → **Settings**), various device settings can be adjusted:

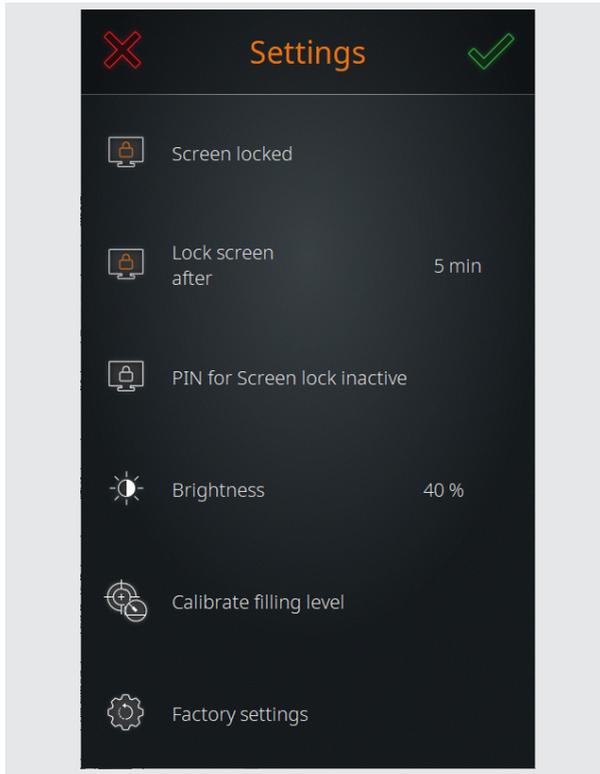The various options are described in detail in the following sections.

Screen lock

The **Screen lock** menu item (**Main menu** → **Settings**) can be used to activate the automatic screen lock for the display (**Screen locked**) or to deactivate it (**Screen unlocked**).

When the screen lock is active, the display locks after a certain period of inactivity (no operator actions on the control panel).

The progress window (**Default**) is displayed in the lock screen. As soon as the display is touched, the unlock slider appears.

The screen lock prevents accidental operation during operation, but does not provide protection against unauthorized access.

To protect the device from unauthorized access, the **PIN for Screen lock** option must be enabled!

Unlock screen

To unlock the screen, touch any point on the screen and swipe the slider from left to right:

Time until automatic screen lock is activated

Using the option **Look screen after . . . min** you can specify after which period of time the automatic screen lock should be activated.

- Touch the function button to open the numeric on-screen keyboard and set the desired time interval in minutes.
- Confirm the entry with .
- To discard the entry and return to the previous menu without making any changes, touch the  button.

Screen lock PIN

The **PIN for Screen lock** option allows you to set a security code that must be entered to unlock the screen. This feature can be used to protect the device from unauthorized access.

- Touch the **PIN for Screen lock** function button to open the on-screen numeric keypad.
- Enter a four-digit PIN and press  to confirm.
- Re-enter the PIN to confirm (security prompt!) and end with  to activate the security code → **PIN for Screen lock active**.
- Touch the **PIN for Screen lock** function button again to disable the security code → **PIN for Screen lock inactive**.

The security code is automatically disabled when the screen lock is disabled.

Display brightness

The brightness of the display can be adjusted via the menu item **Brightness (Main menu → Settings)**.

- Touch the **Brightness** function button to open the on-screen numeric keyboard.
- Enter a percentage value between 10 and 100 (in increments of ten) and press  to confirm.
- To discard the entry and return to the **Settings** menu without making any changes, touch the  button.

Calibrate filling level

The flask level sensor S4 can be calibrated via the menu item **Calibrate filling level (Main menu → Settings)**, see section "Calibrate the flask level sensor" on page 96.

Factory settings

The device can be reset to the factory settings via the menu item **Factory settings (Main menu → Settings)**, see section "Restore factory settings" on page 137.

Select the user language

In the **Languages** menu (**Main menu → Languages**), the user language can be changed at any time. Available languages are English, German, Mandarin, French, Spanish, Portuguese, Italian, Russian and Polish.

- Touch the **Languages** function button.
- The list of all available menu languages opens.
- Touch the entry for the desired language.
- The background color of the selected entry changes to orange.
- Press  to confirm the selection and to return to the start screen of the active mode.
- To discard the entry and to return to the start window of the active mode without making any changes, touch the  button.

Select formats

The following system settings can be made in the **Formats** menu (**Main menu** → **Formats**):

- Set date format (EU or US)
- Set time format (12/24-hour)
- Set number separator (EU or US)
- Set temperature unit (degrees Celsius, Fahrenheit, Kelvin)

- The active selection is highlighted in orange. Touch the desired function button to switch between options. Several changes can be made at the same time!
- Confirm the new settings with . All changes are immediately accepted without a confirmation prompt!
- To discard changes and to return to the previous menu, touch the button.

Automatic drainage of residue

In the **Formats** menu, you can define whether the residue is automatically drained from the evaporation flask (blowing out in overpressure mode), or whether the flask is drained by hand.

A system **with automatic drainage of residue**, in sensor or time mode can be operated automatically in continuous operation (condition: continuous supply with educt!).

Operation with automatic drainage of residue

In operation **with automatic drainage of residue**, the residue is automatically drained under the following conditions:

- Stagnation reached in sensor mode.
- Distillate drain counter expired in sensor mode.
- Refill counter expired in sensor mode.
- Additional evaporation time for distillate drain or refill counter has expired in sensor mode.
- Number of defined steps (filling/evaporation) reached in time mode.

Operation without automatic drainage of residue

In operation **without automatic drainage of residue**, the evaporation flask must be drained by hand.

Selection automatic drainage of residue

- The active selection is highlighted in orange when the **Formats** menu is opened. To change, touch the desired function button **with / without automatic. Residue Extraction**.
- Confirm the new setting with . The change is accepted immediately without a confirmation prompt!
- Touch the button  to discard the change.

Select evaporator type/glassware

In the **Formats** menu, you can define the connected evaporator type including the glassware used. The selection stored in the system appears in the form of function buttons:

- The active selection (in the example: Hei-VAP Industrial with set of glassware R) is highlighted in orange. Touch the desired function button to switch between options.
- Confirm the new setting with . The change is accepted immediately without a confirmation prompt!
- Touch the button  to discard the change.


The selection evaporator type/glassware is necessary for display purposes only, see also section "Control panel" on page 90!

Data Recording

The **Data logging** function can be used to record the individual steps of a process. The function records each individual step with start and end time and stores this information in the Data logging memory. The function must be activated before the process starts:


- To activate data recording, switch from the start/overview page of the active mode to the parameters page (**Default**).
- Touch the **Data logging** button in the footer:

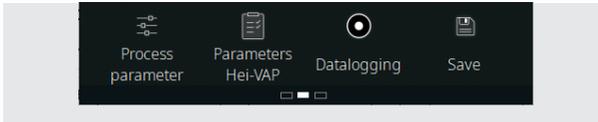- When the Datalogging function is activated, the red recording point appears in the function button and in the progress window display (**default**).

The recorded values can be read out, saved and deleted using the **Data logging** option

- Open the **Main menu** and touch the function button **Data logging**.
- The Data logging list opens:

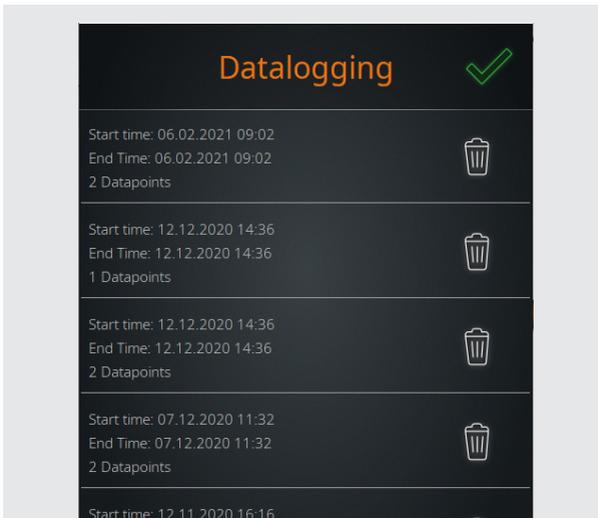- The recorded values show when a process was started and stopped again.

Delete entries

To delete an entry from the list, touch the associated delete button.

Save entries on USB-stick

Recorded entries can be saved on an USB-stick. To do this, connect a VFAT-formatted USB stick to the USB port of the display and open the data list as described.

As soon as a USB stick is connected, a **Save** button (diskette symbol) is assigned to each list entry.

Touch the associated **Save** button to save an entry to the attached USB stick.

PRACTICAL TIP


When recording process data for a long time (several days), large amounts of data accumulate.

It is therefore recommended to save the current status of the data logging memory on a USB stick before starting the process and to remove the existing entries from the internal memory.

System check

All connected sensors and actuators or peripherals are displayed in the **System check** menu. Various tests can be carried out via the **System check** sub-menu to ensure the functional and operational safety of the system:

- Open the **Main menu** and touch the **System check** function button in the footer. The **System check** menu appears:

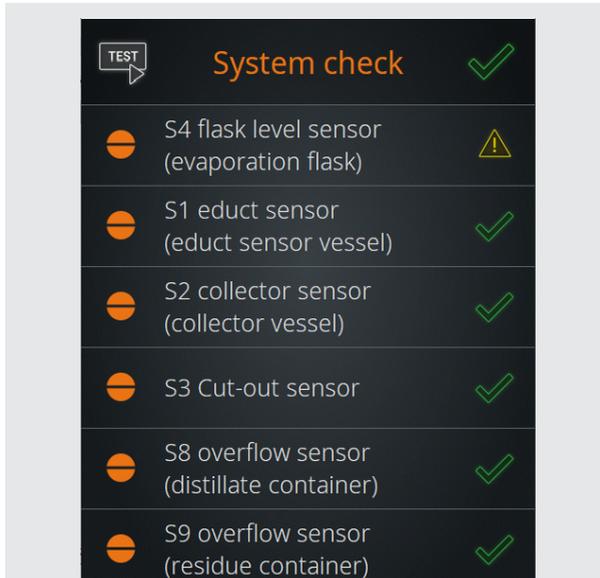- **Error-free** components are indicated by the status display.
- **Faulty** components are marked with , the corresponding entries move to the top of the list (in the example, flask level sensor S4). Verify that the affected component is connected correctly.
- Swipe up or down on the touch panel to scroll through the list.
- Touch the **TEST** function button in the header of the **System check** menu to open the **System check** sub-menu.

CAUTION: Property damage, loss of production

Faulty connections or signaling can lead to malfunctions or uncontrolled switching operations!

Leaks reduce the performance of the entire system.


- Correct any indicated errors and/or detected malfunctions immediately.
- In the event of a continued fault indication/malfunction, contact your local sales representative or our technical service department, see "Contact information Heidolph international" on page 145.

Test valves and signaling

The Test **all actuators** function can be used to test the function of all connected valves and signal lights.

→ Open the **System check** menu as described in section "" on page 109:

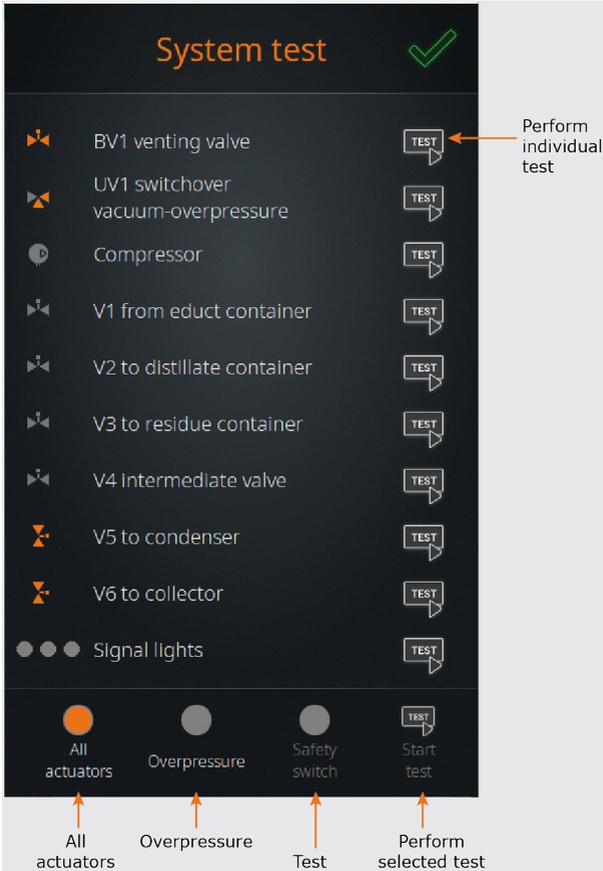Test all valves and signal lights

- In the footer of the **System check** overview, touch the **All actuators** function button to activate the test function (the button icon changes from gray to orange).
- Start the test with **Start test**:
 - During the test, all valves and the signal lights are switched twice in succession (note the color change of the corresponding symbol!). At the end of the test, all switching elements are back in their initial state.
- The test can be stopped at any time by **Stop test**.


Test individual valves and signal lights

- In the **System check** menu, each individual entry is assigned a **TEST** button.
- Touch the **TEST** button of the desired component to start an individual test.
 - The valve or signal light is switched once (note the color change of the corresponding symbol!).
 - Several individual tests can be carried out in parallel.
- Touch the associated **TEST** button again to end the individual test. The individual tests are not automatically terminated!
- Press  to close the menu.

Perform overpressure test

With this test function, the system can be tested for leaks.

- Open the **Check system** menu as described in section "" on page 109.
- Touch the **Overpressure** function button in the footer to activate the test function (button changes to orange).

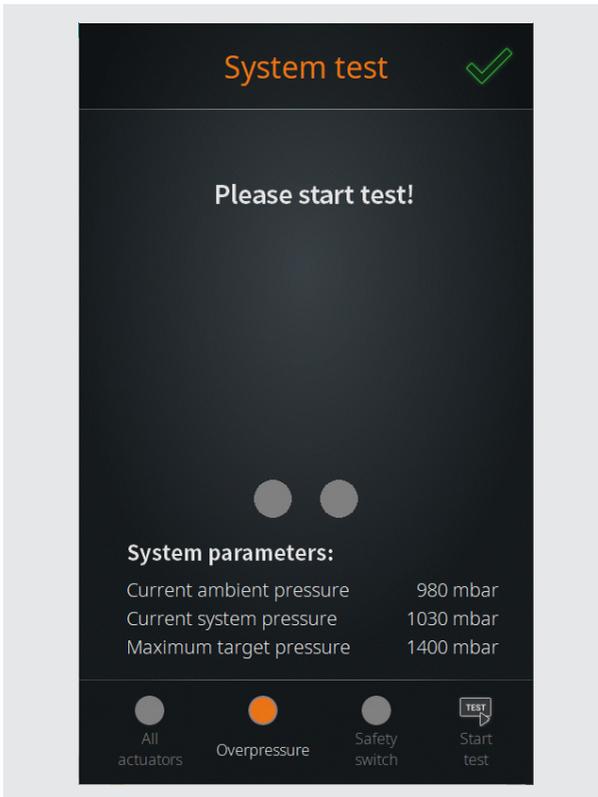WARNING: Risk Of Injury, Damage To Property

Due to faulty components, faulty assembly of components or faulty connections, there is a risk that pressure-applied components will suddenly loosen or rupture during an overpressure test!

- Before performing an overpressure test, ensure that all components are properly installed.
- In particular, make sure that all screw connections are properly tightened.
- Before starting the overpressure test, close all safety covers, hoods and doors of the entire system.
- Ensure that third parties are clear of the danger zone.
- Use the required personal protective equipment (eye protection, hand protection, etc.).

- Touch the **Start test** function button.
- Confirm the subsequent confirmation prompt with  to carry out the overpressure test:
 - During the test, the system pressure is built up to the defined maximum pressure (**Settings** menu). This pressure level must be reached within a maximum of 30 seconds!
 - When the maximum pressure is reached, a measurement is made for a period of 240 seconds to determine whether the system pressure drops below a defined limit value.
 - After this time has elapsed, the system message **Overpressure test successfully completed** is shown on the display, provided that all required values have been achieved.
- If the defined maximum pressure is not reached during the overpressure test or the system pressure falls below a defined limit value, the error message **Overpressure test failed** appears. In this case, follow the instructions on the display.
- The test can be stopped at any time by **Stop test**.
- Press  to close the menu.

Test the safety switch

The function of the integrated safety switch can be tested with the test function **Safety switch**.

The safety switch is activated when a defined system overpressure is reached during compressor operation. When the safety switch is activated, the compressor is shut down so that the system pressure can be released.

- Open the **Check system** menu as described in section "" on page 109.
- In the footer, touch the function button **Safety switch** to activate the test function (button changes to orange).
- Touch the **Start test** function button.
- Confirm the subsequent security query with  to carry out the safety switch test:
 - During the test, the compressor is used to build up pressure in the system. After a certain time has elapsed, the pressure mode valve is automatically switched to vacuum mode by software and the compressor line is blocked abruptly.
- If the safety switch does not trigger, a corresponding error message is displayed and the red point in the middle of the display lights up continuously.
- If the safety switch triggers correctly, the green point in the middle of the display lights up continuously and signals that the safety switch is working properly.

If the safety switch triggers correctly, the standard error message for switching off the system after the safety switch has been triggered also appears!

Follow the instructions on the display.

- The test can be stopped at any time by **Stop test** .
- Press  to close the menu.


Select operating mode

After switching on the device and initializing the device control, the start window of the last active operating mode appears on the display. Touch the **Mode** function button to open the **Applications** menu and to activate the desired mode (see section "Applications menu" on page 102).

Sensor mode

Activate sensor mode

- Open the menu **Applications** and touch the function button **Sensor mode**. The entry is highlighted in orange.
- Confirm the selection with  to activate the mode and open the sensor mode start window:

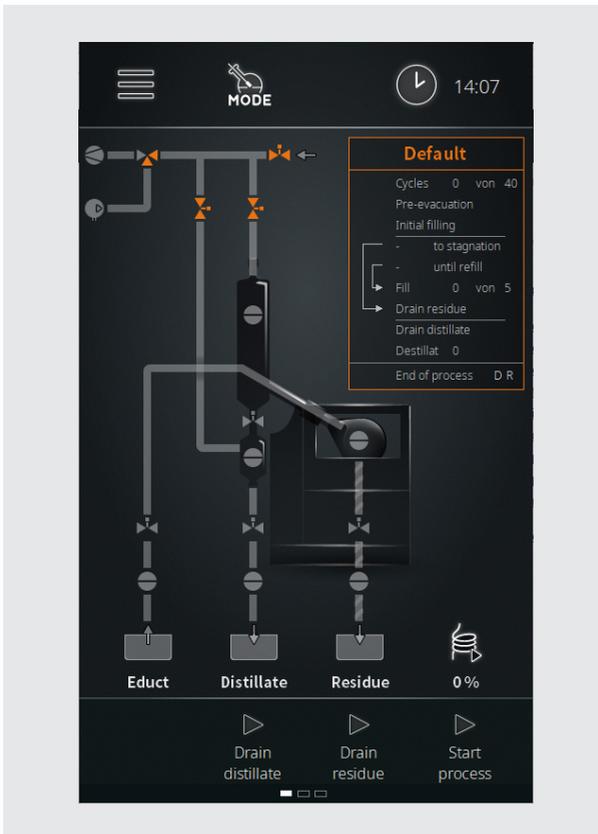In sensor mode, the filling quantity in the evaporation flask is regulated by the flask level sensor S4. The following sensor statuses are possible:

- The filling process is started and runs as long as the sensor does not have contact with the medium → S4 in the status **low**.
- The filling process is stopped when the sensor comes into contact with the medium → S4 in the state **high**.

Prerequisites

The following conditions must be met for operation in sensor mode:

- Heating bath heater of the rotary evaporator is switched on
- Recirculation chiller is switched on
- Vacuum pump is switched on
- Rotary evaporator is under filling vacuum
- Level sensor S4 calibrated
- Educt sensor vessel is filled


Filling the evaporation flask

After the initial filling of the evaporation flask, the filled educt is brought to process temperature.

For the initial filling, it is essential to observe the instructions in section "Initial filling of the evaporation flask" on page 99.

- Touch the **Start process** button in the footer of the start window. The following message appears on the display:

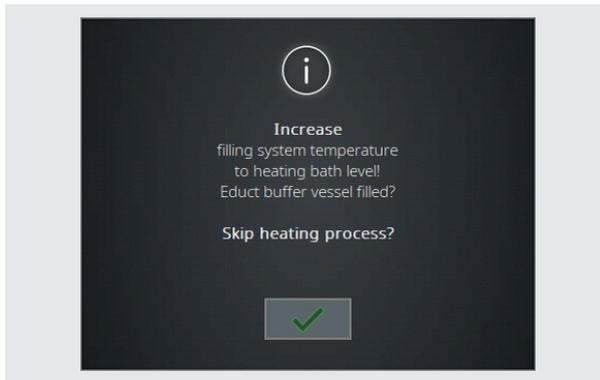
- The filled-in template is saved within a preset time (factory setting: 300 s) to the heating bath temperature.
 - Observe the progress bar in the message box.
 - Use  to skip the heating process.

Note that skipping the heating process will automatically start the evaporation process. This activates the refill, the emergency stop and the feed sensor!

- After the pre-heating time has elapsed or after the step has been skipped, the filling pressure is automatically regulated to the set process pressure. The following message appears on the display:

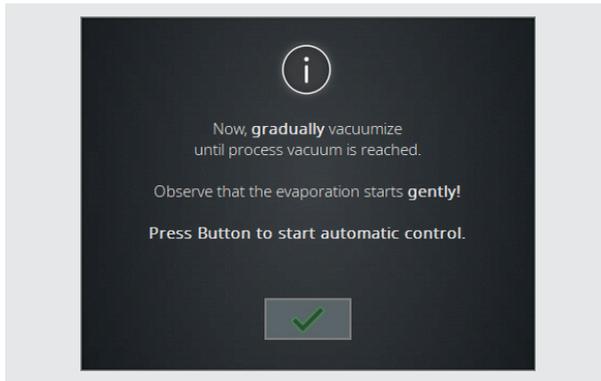- Start the automatic mode with , as soon as the evaporation process starts (condensate is deposited in the glass cooler).

Automatic operation

The process time starts when automatic mode is started.

- As the evaporation starts, the filling level in the evaporation flask decreases:
 - The flask level sensor S4 loses contact with the medium and changes from **high** to **low**. With this status change, the **filling time** is started (→ **low** = start filling time).
- As soon as the defined setpoint for the filling time has been reached (**filling time** expired), the evaporation flask is automatically filled with new educt.
 - The **stagnation time** starts when the level sensor S4 changes from low to high during filling (→ **high**= start stagnation time).
- The filling level in the evaporation flask rises again. As soon as the level sensor S4 detects media contact, the status changes again from **low** to **high**(→ **high** = start stagnation time).
- The continuous evaporation process reduces the filling level in the evaporation flask again:
 - The level sensor S4 loses contact with the medium and changes from **high** to **low** (→ **low** = start filling time).
 - When the status changes again, the **filling time** is started and the **stagnation time** is set to zero.
 - As soon as the defined **filling time** has been reached again, the evaporation flask is filled with new educt.
- The filling level in the evaporation flask rises again and the described process steps are repeated until the end of the process.

The **filling time** defines the time until the next filling of the evaporation flask and thus the amount of educt that is evaporated until the refilling step.

The **stagnation time** runs as soon as the level sensor S4 changes to the high status after filling. If the stagnation time expires completely, the automatic residue emptying is initiated (depending on the set of glassware and the parameterization) or the piston must be emptied manually.

End of process in sensor mode

In sensor mode, the process sequence is stopped under one of the following conditions:

- Educt container runs empty
 - After the **filling time** has elapsed, the evaporation flask is no longer filled → Status change of the feed sensor S1 (**high to low**).
 - Unplanned end is initiated: Additional evaporation time is started (setting via process parameters, see section "Parameter sensor mode" on page 119).
 - After the additional evaporation time has elapsed, the system switches to a safe operating state: Peripheral shutdown (recirculation chiller and vacuum pump off, evaporator in standby).
- Maximum filling rate reached
 - Using the process parameter **max. filling: Additional Evaporation time** a maximum filling rate can be set in order to achieve a defined concentration of the residue (setting via process parameters, see section "Parameter sensor mode" on page 119).
- Maximum number of distillate drainage steps reached
 - Use the process parameter to **empty max. distillate: Additional Evaporation time**, a maximum number of distillate drainages can be specified in order to achieve a defined concentration of the residue (setting via process parameters, see section "Parameter sensor mode" on page 119).
- Residue level in the evaporation flask reached, flask level sensor S4 remains permanently in the status **high** (medium contact), nominal **stagnation time** is reached.
 - In this case, the sensor control is stopped and the residue is drained automatically or manually. The draining process depends on the glassware (with/without automatic drainage of residue) and the default **Drain residue on/off** in the **Default** window (setting via process parameters, see section "Parameter sensor mode" on page 119).
- Maximum filling level in distillate/residue container → Level sensor S8 (distillate container) or S9 (residue container) indicate by a status change from **low to high** that the maximum level has been reached (option).


Both of the level sensors S8 (distillate container) and S9 (residue container) are available as optional accessories and are not included in the standard scope of delivery!

The level sensors prevent overflow of the monitored vessel. When a level sensor is triggered, the system is switched to a safe operating state.

Parameter sensor mode

General parameters

- To adjust the **general parameters**, switch from the start window of the sensor mode to the corresponding **Default** view (tap the default window or swipe from right to left):

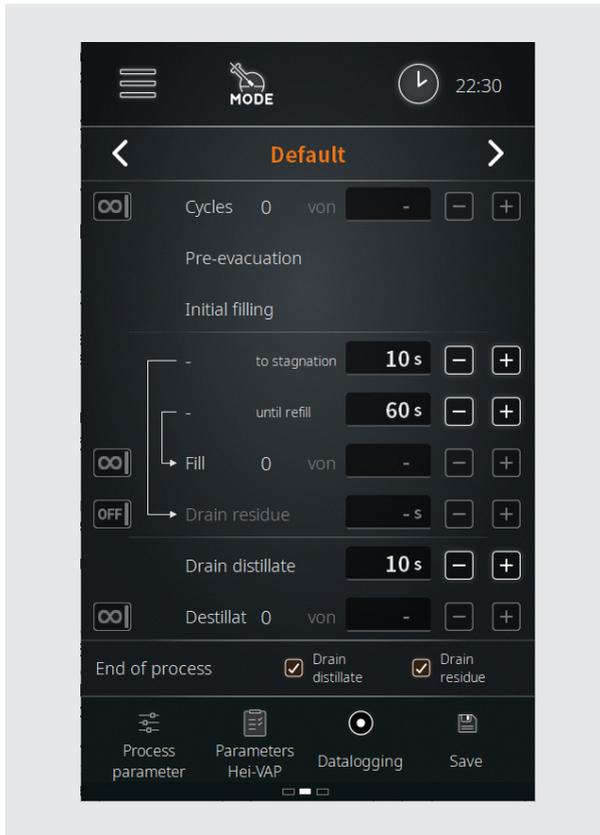- Use the two navigation keys to return to the start window or to continue to the **Recipes** view:

The following parameters can be adjusted in the **Default** view:

| Parameters | Function |
|-------------------------|---|
| Cycles | <p>Maximum number of evaporation cycles</p> <ul style="list-style-type: none"> Selection switch in position ∞ → no limit. Selection switch in position n → define a certain number of cycles in the selection field (adjustment via plus/minus keys). <p>Touch the ∞ / n icon to switch between the two options.</p> |
| Fill | <p>Maximum filling rate</p> <ul style="list-style-type: none"> Selection switch in position ∞ → no limit. Selection switch in position n → define a maximum number of filling operations in the selection field (adjustment via plus/minus keys). <p>Touch the ∞ / n icon to switch between the two options.</p> |
| Drain residue | <p>Residue draining time</p> <ul style="list-style-type: none"> Selection switch in position off → no defined draining time. Selection switch in position on → define a residue draining time in the selection field (adjustment via plus/minus keys). <p>Touch the on / off icon to toggle between the two options.</p> |
| Drain distillate | <p>Distillate draining time</p> <ul style="list-style-type: none"> In the selection field, define a distillate draining time (adjustment via plus/minus keys). |
| Distillate | <p>Maximum number of distillate drainage steps</p> <ul style="list-style-type: none"> Selection switch in position ∞ → no limit. Selection switch in position n → define a maximum number of distillate drainage steps in the selection field (adjustment via plus/minus keys). <p>Touch the ∞ / n icon to switch between the two options.</p> |
| End of process | <p>Specify here whether the distillate and/or residue are drained at the end of the process. Touch the corresponding check box to enable (check mark) or disable (blank) the option.</p> |

The parameters **Fill** and **Distillate** are mutually exclusive! For an evaporation cycle, you can specify either a maximum filling rate **OR** a maximum number of distillate drainage steps.


→ Confirm or cancel all changes with **ok** or **cancel**.

Process parameters

- Touch the **Process parameters** button in the footer to open the menu.
- Tap the button of the desired parameter in the menu and adjust the value using the on-screen keyboard, see also "User input" on page 95.

| Parameters | Function |
|--|---|
| Heating tape | Performance of the heating tape during operation. Range: 10 – 100% in increments of ten. |
| Heating tape automatic | Automatically switch on heating tape yes/no. Range: on / off. |
| Evacuation after draining collector | Evacuation period of the bypass line after draining the collector vessel. Range: 1 – 9999 s. |
| Warm-up time filling in sensor mode | Warm-up time when filling the evaporation flask in sensor mode. Range: 1 – 9999 s. |
| Educt empty: additional Evaporation time | When the educt container runs empty: Additional evaporation time for evaporation of the remaining educt in the evaporation flask. Range: 1 – 9999 s. |

Maximum filling rate: Selectable when a specific number of filling operations is specified in the start window:

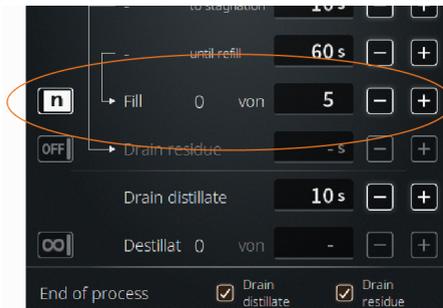
Max. filling: additional Evaporation time

In the example: five fillings, then optionally an additional evaporation time runs. This allows a desired concentration to be achieved!

The parameters **max. filling: Add. Evaporation time** and **max. drain distillate: additional Evaporation time** are mutually exclusive and cannot be activated at the same time!

Selectable if a specific number of distillation drainage steps is specified in the start window:

Max. drain distillate: additional Evaporation time

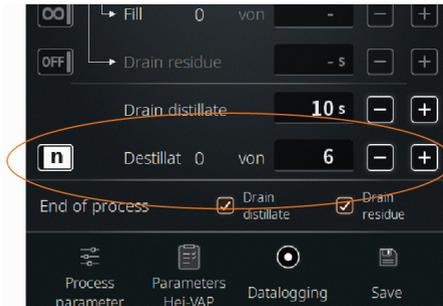In the example: five distillate drainages, then optionally an additional evaporation time runs. This allows a desired concentration to be achieved!

The parameters **max. filling: Add. Evaporation time** and **max. drain distillate: additional Evaporation time** are mutually exclusive and cannot be activated at the same time!

Forced draining collector

Forced draining collector active: Collector will be drained after 1 – 9999 seconds, regardless of the filling level.

Range: 1 – 9999 s.

Activate/deactivate forced draining of the collector (determines previous parameter!).

Range: active / inactive

Standby timer after educt container runs empty

Time until the standby mode is activated after the educt container runs empty.

Range: 1 – 9999 s.

Parameters Hei-VAP

In the sub-menu **Parameters Hei-VAP** various parameters of the connected Hei-VAP can be stored (vacuum, hysteresis, rotational speed, heating bath temperature, recirculation chiller temperature).

These values are for informational purposes only and may be included in a new recipe.

Touch the **Parameters Hei-VAP** button in the footer of the **Default** view to open the sub-menu and enter the desired values of the connected rotary evaporator by hand.


Data logging

See section "Data Recording" on page 108.

Recipes

By touching the **Save** button, the actual settings and values can be stored in the system as a recipe.

→ Touch the **Save** button to open the on-screen keyboard (see also section “User input” on page 95):

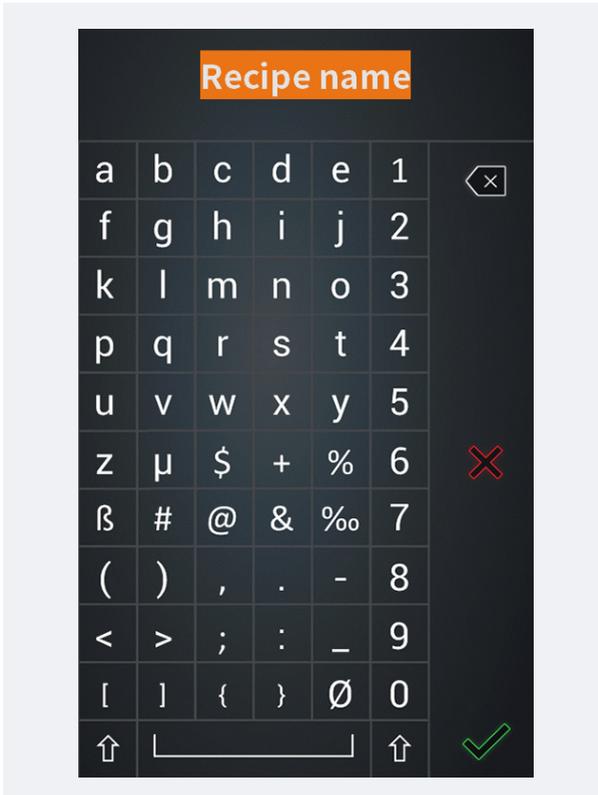→ Use the on-screen keyboard to enter any recipe name.

It is not possible to assign recipe names twice! If a desired recipe name is already stored in the system, a corresponding warning message appears.

In this case, you can enter a different name or update the existing recipe with the new values.

→ Confirm your entry with  to store the new recipe or the new settings in the recipe list.

Recipe list

To open the recipe list, switch from the start window of the sensor mode to the **Recipes** view (two swipes from right to left, or use the navigation button from the **Default** view).

Here, new recipes can be created (**New** button) and existing recipes can be loaded (**Load** button) or deleted (**Delete** button).


The entire recipe list can be saved on a USB stick. Recipe lists can also be loaded from a USB stick to the device memory.

To do this, connect a USB stick to the USB port on the control panel and follow the instructions on the display.

Time mode

Activate time mode

- Open the menu **Applications** and touch the function button **Time mode**. The entry is highlighted in orange.
- Confirm the selection with  to activate the mode and open the time mode start window:

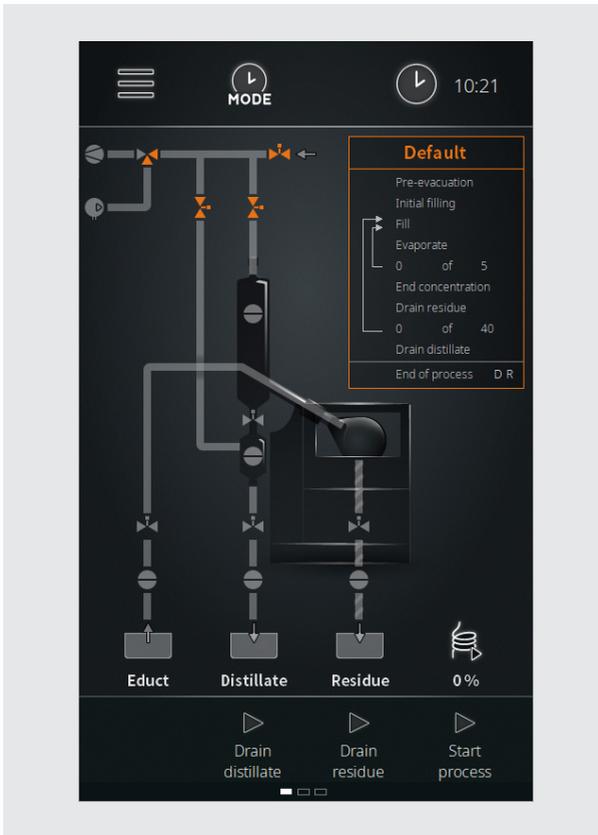In the time mode, the filling of the evaporation flask, the distillation and the draining of the distillate and residue are carried out in a purely time-controlled manner:

- Within a cycle, a desired amount of concentrate is built up as a residue in the evaporation flask in the time mode.
- A cycle can consist of 1 – 9999 steps.
- One step consists of the **filling** and **evaporation** processes.
- In each step, the desired end concentration is reached.

Prerequisites

The following conditions must be met for operation in sensor mode:

- Heating bath heater of the rotary evaporator is switched on
- Recirculation chiller is switched on
- Vacuum pump is switched on
- Rotary evaporator is under filling vacuum
- Process parameters are set correctly on the rotary evaporator
- Level sensor S4 calibrated

Before the initial filling of the evaporation flask or before the first process step, it must be ensured in time mode that all process parameters (vacuum, temperature, speed, cooling) are set correctly on the rotary evaporator.

Otherwise, there may be inaccuracies in the recipe processing!

Procedure

- Initial filling in time mode: This step is only necessary to prevent drying or crystallization of the educt. The activation is carried out via the parameter **Initial filling** in the **Default** view, see the following section "Time mode parameters" on page 127. For the initial filling, it is essential to observe the instructions in section "Initial filling of the evaporation flask" on page 99.
- In time mode, a defined quantity is filled into the evaporation flask and immediately vaporized within a specified time:
 - The filling quantity is defined by the **filling time**.
 - The **evaporation time** is defined by the user.
- Due to the **pre-evacuation** of the evaporation flask, the evaporation process starts immediately while filling. After processing a step, the evaporation flask is filled again and the following step is processed:
 - The number of steps in a cycle is parameterized by the user. This value depends on the desired target quantity, which is to form in layers in the evaporation flask.
- After processing the specified number of steps, a cycle has been run through and the residue is pressed out of the evaporation flask by means of overpressure.
 - The number of cycles is parameterized by the user.
- At the end of the last cycle, all connected devices are switched off (operating safe state, end of process).


In time mode, the system containers are monitored by sensors! The end process can be initiated and/or an ongoing process can be terminated and the evaporator can be put into a safe operating state when a corresponding signal is emitted by a sensor.


Time mode parameters

General parameters

- To adjust the **general parameters**, switch from the start window of the time mode to the corresponding **Default** view (tap the default window or swipe from right to left):

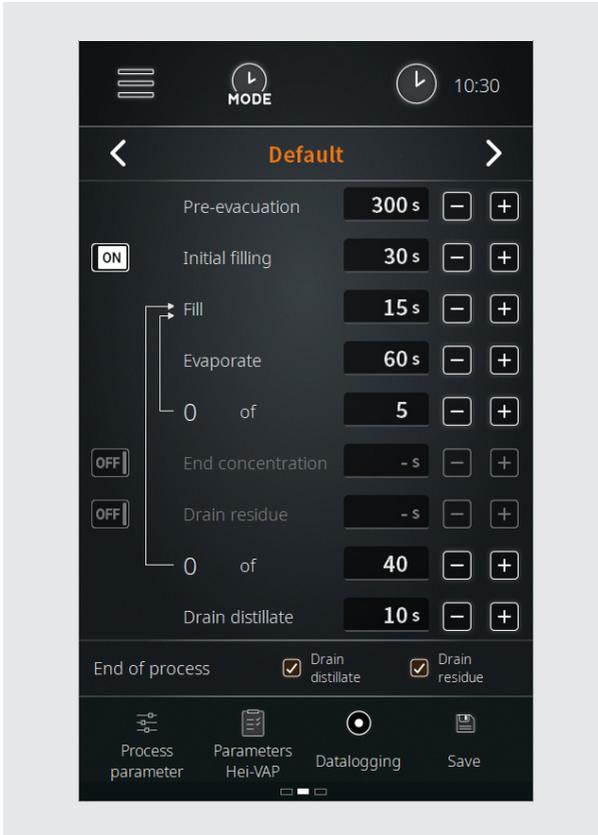- Use the two navigation keys to return to the start window or to continue to the **Recipes** view:

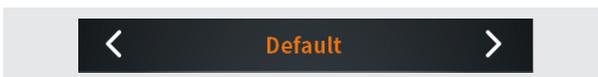The following parameters can be adjusted in the **Default** view:

| Parameters | Function |
|--------------------------|--|
| Pre-evacuation | <p>Adjustable time for the pre-evacuation phase Before starting the process loops evaporation – filling – draining, the system is pre-evacuated for the time set here.</p> |
| Initial filling | <p>Adjustable time for initial filling Before starting the process loops evaporation – filling – draining, the system can be initially filled. An initial filling ensures that a minimum quantity is filled before starting the evaporation, avoiding thus the recrystallization of the educt.</p> <ul style="list-style-type: none"> ▪ Selection switch in position off → no initial filling. ▪ Selection switch in position on → specify a time period for the initial filling in the selection field (adjustment via plus/minus keys). <p>Touch the on / off icon to toggle between the two options.</p> |
| Fill | <p>Filling time setting Adjust the filling time as required using the plus/minus keys.</p> |
| Evaporate | <p>Evaporation time setting Adjust the evaporation time as required using the plus/minus buttons.</p> |
| End concentration | <p>Additional evaporation step to reach a specific end concentration Activate this option to pass an additional evaporation step after passing through the process loops to achieve a desired end concentration or to remove any residual moisture from the residue in solid-liquid mixtures.</p> <ul style="list-style-type: none"> ▪ Selection switch in position off → no additional evaporation step. ▪ Selection switch in position on → specify a time period for the one additional evaporation step in the selection field (adjustment via plus/minus keys). <p>Touch the on / off icon to toggle between the two options.</p> |
| Drain residue | <p>Residue draining time</p> <ul style="list-style-type: none"> ▪ Selection switch in position off → no defined draining time. ▪ Selection switch in position on → define a residue draining time in the selection field (adjustment via plus/minus keys). <p>Touch the on / off icon to toggle between the two options.</p> |
| Drain distillate | <p>Distillate draining time</p> <ul style="list-style-type: none"> ▪ In the selection field, define a distillate draining time (adjustment via plus/minus keys). |

End of process

Specify here whether the distillate and/or residue are drained at the end of the process. Touch the corresponding check box to enable (check mark) or disable (blank) the option.

→ Confirm or cancel all changes with **Accept** or **Cancel**.

Process parameters

- Tap the **Process parameters** button in the footer.
 - The process parameters menu opens.
- Touch the function button of the desired parameter and adjust the value using the on-screen keyboard, see also "User input" on page 95.

| Parameters | Function |
|---|--|
| Heating tape | Performance of the heating tape during operation. Range: 10 – 100% in increments of ten. |
| Heating tape automatic | Automatically switch on heating tape yes/no. Range: on / off. |
| Evacuation after draining collector | Evacuation period of the bypass line after draining the collector vessel. Range: 1 – 9999 s. |
| Educt empty: additional Evaporation time | When the educt container runs empty: Additional evaporation time for evaporation of the remaining educt in the evaporation flask. Range: 1 – 9999 s. |
| Forced draining collector | Forced draining collector active: Collector will be drained after 1 – 9999 seconds, regardless of the filling level. Range: 1 – 9999 s. Activate/deactivate forced draining of the collector (determines previous parameter!). Range: active / inactive |
| Standby timer after educt container runs empty | Time until the standby mode is activated after the educt container runs empty. Range: 1 – 9999 s. |


- **Parameters Hei-VAP**, see section "Parameters Hei-VAP" on page 122.
- **Datalogging**, see section "Data Recording" on page 108.
- **Save**, see section "Recipes" on page 123.

Rinse mode

Activate the rinse mode

- Open the menu **Applications** and touch the function button **Rinse mode**. The entry is highlighted in orange.
- Confirm the selection with  to activate the mode and open the rinse mode start window:

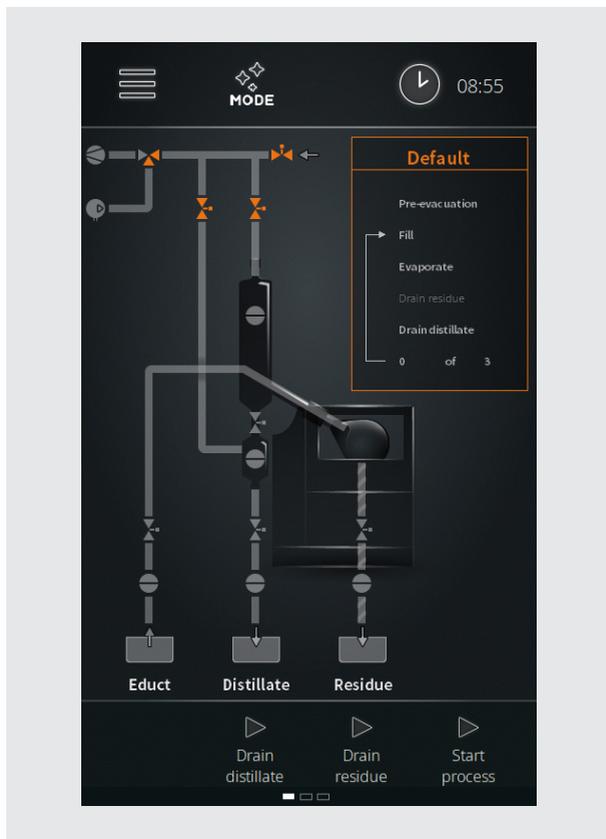The rinse mode allows the fill and drain lines to be flushed with a suitable solvent. The flushing process removes residues from previous processes and impurities. When rotation is activated, the evaporation flask is also cleaned.

Depending on the degree of contamination, several cycles should be carried out!

To clean the complete system including the distillation area and the condenser, an evaporation process with a suitable solvent is necessary. Adjusting the parameters appropriately, the distillate produced this way will moisten all glass parts, thereby cleaning them.

PRACTICAL TIP

- Select a vacuum value adequate for the solvent and set the hysteresis to 50 mbar.
- A large pressure range is passed through by the large hysteresis.
- All glass parts are first wetted, after which drops are formed, which then run off.
- As a result, the glass parts in the distillate area are cleaned within a short time.

Prerequisites

The following conditions must be met for rinsing the distillate tube:

- Heating bath of the rotary evaporator is switched on
- Recirculation chiller is switched on
- Vacuum pump is switched on


Procedure

- Depending on the selection, the various lines of the system are flushed.
- To fill the evaporation flask, see also section "Commissioning" on page 96.


Rinse mode parameters

General parameters

- To adjust the **general parameters**, switch from the start window of the rinse mode to the corresponding **Default** view (tap the default window or swipe from right to left):

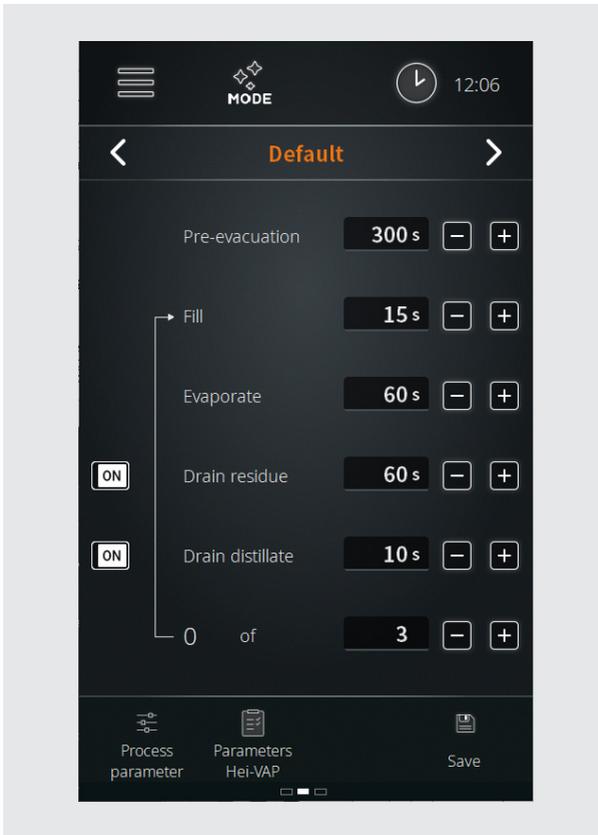- Use the two navigation keys to return to the start window or to continue to the **Recipes** view:

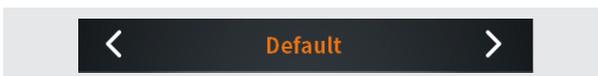The following parameters can be adjusted in the **Default** view:

| Parameters | Function |
|-------------------------|--|
| Pre-evacuation | <p>Adjustable time for the pre-evacuation phase Before starting the process loops evaporation – filling – draining, the system is pre-evacuated for the time set here.</p> |
| Fill | <p>Filling time setting Adjust the filling time as required using the plus/minus keys.</p> |
| Evaporate | <p>Evaporation time setting Adjust the evaporation time as required using the plus/minus buttons.</p> |
| Drain residue | <p>Residue draining time</p> <ul style="list-style-type: none"> ▪ Selection switch in position off → no defined draining time. ▪ Selection switch in position on → define a residue container draining time in the selection field (adjustment via plus/minus keys). <p>Touch the on / off icon to toggle between the two options.</p> |
| Drain distillate | <p>Distillate draining time</p> <ul style="list-style-type: none"> ▪ Selection switch in position off → no defined draining time. ▪ Selection switch in position on → define a distillate container draining time in the selection field (adjustment via plus/minus keys). <p>Touch the on / off icon to toggle between the two options.</p> |

Process parameters

Tap the **Process parameters** button in the footer.

- The process parameters menu opens.

→ Touch the function button of the desired parameter and adjust the value using the screen keyboard, see also "User input" on page 95.

| Parameters | Function |
|--|--|
| Heating tape | Performance of the heating tape during operation. Range: 10 – 100% in increments of ten. |
| Heating tape automatic | Automatically switch on heating tape yes/no. Range: on / off. |
| Standby timer after educt container runs empty | Time until the standby mode is activated after the educt container runs empty. Range: 1 – 9999 s. |

- **Parameters Hei-VAP**, see section "Parameters Hei-VAP" on page 122.
- **Save**, see section "Recipes" on page 123.

Manual mode

Activate manual mode

- Open the menu **Applications** and touch the function button **Manual mode**. The entry is highlighted in orange.
- Confirm the selection with  to activate the mode and open the manual mode start window:


In manual mode, single batches can be processed. For automatic filling of the evaporation flask and draining of the distillate, the **Drain distillate** and **Start process** buttons are available on the operating panel.

Operating parameters must be individually defined on the connected rotary evaporator.

Procedure

- Touch the button **Start process**.
 - The vacuum pump is switched on, the pressure in the evaporation flask is lowered to process pressure.
 - In the footer, the **Fill educt** and **Stop process** buttons are displayed

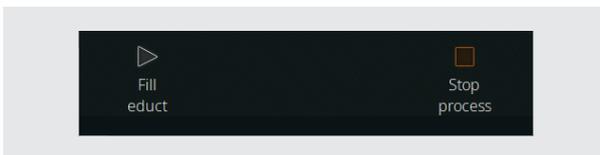- Touch the button **Fill educt** to start the automatic filling of the evaporation flask.

By briefly tapping the **Fill educt** button, a small amount can be injected into the evaporation flask.

Touch the **Fill educt** button for approximately 2 seconds (observe button progress) to fill the evaporation flask continuously. The filling has to be stopped manually!

- Touch the button **Fill educt** again to stop filling.
- Touch the **Stop process** button to stop the process (end process).

When the process is ended, the **Drain distillate** button appears on the control panel.

- Touch the button **Drain distillate** → The distillate is drained (default: approx. 30 s).
- Touch the key **Drain distillate** → Drainage of the distillate is stopped.

The evaporation flask must be drained by hand.

Please observe the relevant information and instructions in the operating instructions for the rotary evaporator!

View error list

- To view the error list, open the **Main menu** and touch the **Error List** function button in the footer.
- To view an individual error message, tap the corresponding entry:
 - The error message appears in plain text on the display.
 - The entries **cannot** be deleted for service reasons!
- Close the error list with . The display returns to the **Main menu**.

Restore factory settings

Use the following procedure to restore the factory settings:

- From the **Main menu**, open the **Settings** menu.
- Touch the **Factory settings** function button.
- The following prompt appears on the display.

- Confirm the security prompt with  to restore the factory settings. The device must be restarted!
- To cancel the reset and return to the previous menu without making any changes, touch the button .

Access system information


The **Information** view displays all version and license information for your automatic module.

- To view system information, open the **Main menu** and touch the **Info** function button in the footer.
- Use the **Versions** and **Licenses** buttons to access the desired information about the software and firmware for your product.
- Scan the stored QR code with a suitable device to access additional product information and relevant documents via our homepage.
- Close the view with . The display returns to the **Main menu**.


Firmware and software update

The system is updated using an update file stored on a USB stick. The upgrade includes the system's software and firmware by default.

- Open the safety cover of the USB port on the control panel and insert the USB stick with the update file.
- Open the **Main menu**:

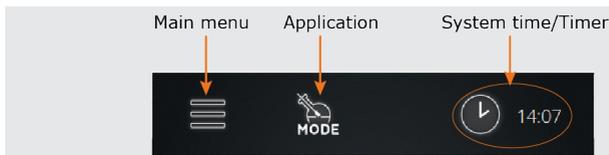- Tap the **System update** button at the bottom of the main menu (this button appears only if a USB stick with a valid update file is detected on the USB port!):

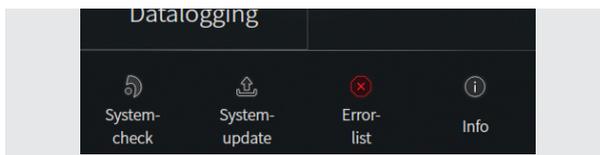- Press  to confirm the security prompt and to update the device software.
- Follow the instructions on the display
- To cancel the operation and return to the previous menu without making any changes, touch the button .
- Restart the device after the software update.
- After the restart, the firmware of the device can be updated. To do this, also confirm the corresponding security prompt with .
- Watch the progress bar and restart the device after the firmware update!
- After the update, unplug the USB stick and put the security cover back into the USB port.

Troubleshooting

| Error message | Possible remedy |
|--|--|
| Cable breakage or sensor defect | Check the cable of the displayed component. |
| Short circuit | Check the cable connections of the displayed component. |
| Values out of range | Correct the values of the displayed parameter. |
| Compressor has been stopped | Excessive pressure from compressor, reduce pressure. |
| Pressure sensor failed! | Check pressure sensor for function, system check. |
| Error due to barometer! | Check the function of the barometer. |
| Temperature sensor failed! | Check temperature sensor for function. |
| Emergency stop has been triggered! | Eliminate the cause of the emergency stop from the external emergency stop switch. |
| Educt container S1 is empty! | Refill educt. |
| Overflow sensor S8 to the distillate container has triggered! | Drain the distillate container. |
| Overflow sensor S9 to the residue container has triggered! | Drain the residue container. |
| S3 cut out sensor has triggered! | Eliminate the cause of the emergency stop on the rotary evaporator. |
| Failure message from chiller | Observe the indicator on the recirculation chiller. |
| Invalid sensor plugged in! | Check the coding of the sensor, check that the sensor is in perfect condition. |
| Invalid valve plugged in! | Check the coding of the valve, check that the valve is in good condition. |
| Communication error! | Check the cable connection and connectors of the displayed components. |
| Operation mode cannot be started! | Open the error list, remove the indicated error cause(s). |

| | |
|---------------------------------------|--|
| Sensor mode cannot be started! | Check the flask level sensor S4 for correct connection. Check filling level in the evaporation flask. Flask level sensor not calibrated? |
| ... has been plugged in twice. | Identical sensors or actuators plugged in, remove one unit. |
| Set switchpoint failed! | Switching error, contact technical service! |
| Calibration failed! | Repeat calibration procedure. |
| Update failed! | Check USB stick for correct connection, restart update. |

All error messages and warning messages are shown in plain text on the display.

Follow the instructions on the display.

In case of recurring errors, please contact the responsible sales department or our technical service. Contact address see "Contact information Heidolph international" on page 145.

Technical Specifications

| General device data | |
|---------------------------------------|---|
| Model | Automatic module Hei-VOLUME Distimatic Pro |
| Dimensions (W× H× D) | 186 × 429 × 521 mm, including control panel |
| Weight | approx. 17 kg |
| IP code housing | IP42 |
| IP code control panel | IP42 |
| Display | Touch panel, 7 Inch |
| Electrical data | |
| Rated voltage | 100 – 240 V AC; 50/60 Hz |
| Connection | L+N+PE |
| Protection class | I |
| Overvoltage category | II |
| Degree of pollution | 2 |
| Power input | 1,500 W max. at 100 V. 3,450 W max. at 230 V. |
| Fuse | 15 A |
| Standard Plug-In Cards | |
| Valve/Sensor IO | Coded connection sockets for valves and sensors |
| Power Control | 1×mains voltage input (plug-in connector C14) |
| | 1×mains voltage output (cold appliance socket F) |
| Permissible ambient conditions | |
| Operating temperature | 5 °C – 31 °C up to 80 % rel. humidity |
| | 32 °C – 40 °C up to 50 % rel. humidity |
| Maximum height above sea level | 2,000 m |
| Storage temperature | 5 °C – 31 °C up to 80 % rel. humidity |
| | 32 °C – 40 °C up to 50 % rel. humidity |

Scope of delivery

| Item | Quantity | Product no. |
|--|----------|----------------|
| Automatic module Distimatic Pro 24/7 Industrial, including control box and set of glassware, with automatic residue extraction | 1 | 591-61112-00-0 |
| Automatic module Distimatic Pro Industrial, including control box and set of glassware, without automatic residue extraction | 1 | 591-61212-00-0 |
| Operating instructions DE/EN | 1 | 01-005-006-51 |
| Guarantee registration / Declaration of no objection | 1 | 01-006-002-78 |
| EC Declaration of Conformity | 1 | 01-001-025-06 |

Further information, in particular about the available expansion modules and the available accessories, can be found on our website at www.heidolph.com!

Device service

When carrying out service work on the device (cleaning, maintenance, repair), observe the general instructions and safety information described in this section.

WARNING: Danger of electric shock

Live components are installed inside the device.

When opening the device, there is a risk of touching live components.

- Switch the device's main switch off and disconnect it from the power supply before carrying out maintenance work, cleaning, or repairs.

Penetrating liquid poses the danger of an electric shock.

- When cleaning, avoid the penetration of liquids.

General cleaning instructions

Wipe all surfaces and the control panel with a damp cloth if necessary. Persistent contamination can be removed with mild soapy water.

CAUTION: Damage to the device

Improper cleaning can damage the surfaces of the device.

Penetrating liquid can damage the electronic components inside the device.

- Clean the device's surfaces with a soft, lint-free and only slightly moistened cloth.
- Never use any aggressive or abrasive cleaning agents or aids.

Repairs

Repairs to the device may only be carried out by authorized skilled experts!

Unauthorized repairs during the warranty period will result in the loss of the warranty claim.

The owner is solely liable for damage caused by unauthorized repairs.

In case of repair contact an authorized dealer or our technical service, see "Contact information Heidolph international" on page 145.

Include the completed declaration of no objection with every device return, see "Declaration of no objection" on page 146.

Maintenance

The device contains no user-serviceable components. If necessary, in the event of abnormal operating behavior such as excessive noise or heat generation, for example, contact our technical service, see "Contact information Heidolph international" on page 145.

Dismantling

Observe the instructions contained in the assembly instructions supplied when dismantling.

Disposal

- When disposing of the device, observe the provisions of the WEEE Directive 2012/19/EU and its transposition into national law in the country of use.
 - When disposing of portable batteries, observe the provisions of the European Battery Directive 2013/56/EU and their transposition into national law in the country of use.
 - Check the device and all components for residues of substances that are hazardous to health, the environment, and biohazardous before disposing.
 - Properly remove and dispose residues of substances that are hazardous to health, the environment and biohazardous!
-

Contact information Heidolph international

Heidolph Instruments North America

Phone: 1-866-650-9604
E-mail: service@heidolph.com
www.heidolphNA.com

Heidolph Instruments United Kingdom

Phone: 01799 - 5133-20
E-mail: service@radleys.co.uk
www.heidolph-instruments.co.uk

Local distributors

To find your local distributor please visit www.heidolph.com

Warranty Statement

Heidolph Instruments assumes a three-year warranty against material and manufacturing defects.

Excluded from the warranty are glass and wear parts, transport damage, and damage resulting from improper handling or non-intended use of the product.

For registered products, the warranty period begins on the date of purchase. Register the product with the enclosed warranty card or on our homepage www.heidolph.com.

For non-registered products, the warranty period begins with the date of the serial production (to be determined by the serial number).

In the event of material or manufacturing defects, the product will either be repaired or replaced free of charge within the warranty period.

Declaration of no objection

Enclose the declaration of no objection, duly completed, with your device return.
Submissions without a declaration of no objection cannot be processed!

1. Information on the device

Item number: _____
 Serial number: _____
 Reason for submission: _____

2. Does this product pose a risk to people and/or the environment due to its use for processing substances that are hazardous to health, the environment and/or biohazardous? Please mark with a cross and complete the information!

| | | |
|-----------|------------|--|
| | | If YES , with which substances did the device come into contact? |
| | | _____ |
| | | _____ |
| NO | YES | |
| | | If NO , what cleaning and/or decontamination measures were carried out? |
| | | _____ |
| | | _____ |

3. Information on the client/sender:

Name, first name: _____
 Company/institution: _____
 Department/working group: _____
 Address: _____
 Zip code, city: _____
 Country: _____
 Phone: _____
 E-Mail: _____

4. Legally binding declaration

With his signature, the client/submitter declares the completeness and correctness of his details. Missing or incorrect information obligates to compensation for damages.

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With this declaration, we confirm (according to current knowledge) that all sold laboratory devices by Heidolph Instruments GmbH & Co. KG fulfill the requirements of the EU directive 2011 / 65 / EU (RoHS) and the amended of directive 2015 / 863.

All devices are compatible with the requirement of the RoHS for the following materials:

Max. 0,01% of the weight in homogeneous material for cadmium and max. 0,1 % of the weight in homogeneous material for lead, mercury, hexavalent chromium, polybrominated biphenyl, polybrominated diphenyl ether, Di (2-ethylhexy l) phthalate, butyl benzyl phthalate, dibutyl phthalate, diisobutyl phthalate.

In the case of individual assemblies, maximum concentrations maybe exceeded within the permissible exceptions of the Directive.

Schwabach, 22.07.2019

Stefan Peters
Vice President Marketing, Innovation & Technology

Marcell Sarré
Vice President Quality Management & Technical Service

Zchnng.-Nr.

China RoHS Certification

China RoHS DECLARATION OF CONFORMITY

Heidolph Instruments GmbH & Co.KG has made reasonable efforts to ensure that hazardous materials and substances may not be used in its products.

In order to determine the concentration of hazardous substances in all homogeneous materials of the subassemblies, a "Product Conformity Assessment" (PCA) procedure was performed. As defined in GB/T 26572 the "Maximum Concentration Value" limits (MCV) apply to these restricted substances:

- Lead (Pb): 0.1%
- Mercury (Hg): 0.1%
- Cadmium (Cd): 0.01%
- Hexavalent chromium (Cr(VI)): 0.1%
- Polybrominated biphenyls (PBB): 0.1%
- Polybrominated diphenyl ether (PBDE): 0.1%

Environmental Friendly Use Period (EFUP)

EFUP defines the period in years during which the hazardous substances contained in electrical and electronic products will not leak or mutate under normal operating conditions. During normal use by the user such electrical and electronic products will not result in serious environmental pollution, cause serious bodily injury or damage to the user's assets.

The Environmental Friendly Use Period for Heidolph Instruments GmbH & Co.KG products is 25 years.

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This table is created according to SJ/T 11364-2014

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Note: Table applies to all products. Some of the components or parts listed above may not be part of the enclosed product.

- O: 表示该有毒有害物质在该部件所有均质材料中的含量均在GB/T 26572规定的限量要求以下。
- O: Indicates that the above mentioned hazardous substance contained in all homogeneous materials of the part is below the required limit as defined in GB/T 26572.
- X: 表示该有毒有害物质至少在该部件某一均质材料中的含量超出GB/T 26572规定的限量要求。
- X: Indicates that the above mentioned hazardous substance contained in at least one of the homogeneous materials of this part is above the required limit as defined in GB/T 26572.

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Apart from the disclosures in the above table, the subassemblies are not intentionally manufactured or formulated with lead (Pb), mercury (Hg), cadmium (Cd), hexavalent chromium (CrVI), polybrominated biphenyls (PBB), and polybrominated diphenyl ethers (PBDE).

Products manufactured by Heidolph Instruments GmbH & Co.KG may enter into further devices or can be used together with other appliances.

With these products and appliances in particular, Heidolph Instruments GmbH & Co.KG will not take responsibility for the EFUP of those products and appliances.

Place, date Schwabach, 19.09. 2019

Wolfgang Jaenicke
Chief Executive Officer CEO

Marcell Sarré
Vice President
Quality Management &
Technical Service

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Subject to change without notice. The printed version of this operating instructions is not regularly updated. The latest issue of this document can be found by visiting our homepage.

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© Heidolph Instruments GmbH & Co. KG